



**Iridex Pascal®
(532 nm / 577 nm)**

**Ophthalmic Scanning Laser System
Operator Manual**

**Iridex Corporation
1212 Terra Bella Avenue
Mountain View, CA 94043 USA
Office: +1.650.940.4700
Fax: +1.650.962.0486
www.iridexretina.com**

Caution: "Federal law restricts this device to sale by or on the order of a physician"

This manual is copyrighted with all rights reserved. Under the copyright laws, this manual may not be copied in whole or in part or reproduced in any other media without the express written permission of Iridex Corporation. Permitted copies must carry the same proprietary and copyright notices as were affixed to the original under the law. Copying includes translation into another language.

Please note that while every effort has been made to ensure that the data given in this document is accurate, the information, figures, illustrations, tables, specifications and schematics contained herein are subject to change without notice.

Iridex[®], the Iridex - Logo[™] and PASCAL[®] are trademarks or registered trademarks of Iridex Corporation.

Iridex Pascal is trademarks of Iridex Corporation.

Iridex Pascal refers to a family of laser products: single wavelength 532 nm and single wavelength 577 nm.

Iridex Pascal Laser Systems Operator Manual DC-04076 Rev F 05.29.2026 © 2026 by Iridex Corporation. All rights reserved.

Introduction

Intended Use

Iridex Pascal is a laser system console with an integrated slit lamp. The system connects to the slit lamp to enable laser energy to be delivered through the slit lamp illumination path. The system may be used for standard single shot photocoagulation and laser scanning patterns.

The system enables the physician to deliver multiple laser spots with a single footswitch depression by automating the emission of laser light. The aiming beam displays the pattern, allowing the physician to place it in the appropriate location.

Indications For Use

Iridex Pascal is intended for use to perform single-spot photocoagulation in the posterior segment (retina, choroid) and in the anterior segment (iris, trabecular meshwork) as well as pattern-scanning photocoagulation in the non-macular retina of the eye. Single-spot delivery may be performed using a slit lamp biomicroscope or an indirect ophthalmoscope. Pattern delivery may be performed using a slit lamp biomicroscope.

The system is intended for use by trained Ophthalmologists for diagnosis and treatment of ocular pathology in both the posterior and anterior segments. Intended for use in the posterior segment to perform retinal photocoagulation, panretinal photocoagulation, focal photocoagulation and grid photocoagulation for vascular and structural abnormalities of the retina and choroid including:

(532 nm)

- proliferative and non-proliferative diabetic retinopathy
- macular edema
- choroidal neovascularization associated with wet age-related macular degeneration
- age-related macular degeneration
- lattice degeneration
- retinal tears and detachments

(577 nm)

- proliferative and non-proliferative diabetic retinopathy
- macular edema
- choroidal neovascularization associated with wet age-related macular degeneration
- age-related macular degeneration
- lattice degeneration
- retinal tears and detachments

Intended for use in the treatment of ocular pathology in the anterior segment including: (532 nm and 577 nm)

- iridotomy
- trabeculoplasty

Purpose Of This Manual

This manual outlines Iridex Pascal, including operating procedures, troubleshooting, maintenance and cleaning. Please refer to this operator manual for instructions and guidance on how to appropriately use your laser.

Before using the system, carefully read the “[CAUTIONS & WARNINGS](#)” and “[General Safety and Regulatory Information](#)” sections to familiarize yourself with the operation of the system.

Attention

Iridex accepts full responsibility for safety, reliability and performance of the device only if:

- Service, readjustments, modifications and/or repairs are performed exclusively by Iridex Corporation-certified personnel.
- The electrical installation of the treatment room complies with the applicable IEC, CEC and NEC requirements.

The warranty is void if any of these warnings are disregarded.

Iridex Corporation reserves the right to make changes to the device(s) herein. Device(s), therefore, may not agree in detail with the published design or specifications. All specifications are subject to change without notice.

For questions regarding your laser, please contact Iridex Corporation or your local Iridex Corporation representative.

Iridex Corporation Contact Information



Iridex Corporation
1212 Terra Bella Avenue
Mountain View, California 94043-1824 USA
Telephone: (650) 940-4700
(800) 388-4747 (US only)
Fax: (650) 962-0486
Technical Support: (650) 962-8100
techsupport@Iridex.com

Table of Contents

Introduction.....	3
Intended Use.....	3
Indications For Use	3
Purpose Of This Manual	4
Attention.....	4
Cautions and Warnings.....	11
Disclaimer	11
Operator Manual Symbol Definitions	11
General Safety and Regulatory Information.....	11
Ocular Protection	12
Laser Safety Eyewear.....	13
Electrical Hazards	14
Non-Homogeneous Pattern Delivery	15
Fiber Optic Cable Assembly	15
Pattern Titration	15
Titration	16
Fire Hazard	16
Protecting Non-target Tissues	16
Operational Safety	17
Additional Safety Considerations	18
Safety Compliance	19
Regulatory Compliance Safety Features	19
Key Lock Switch.....	19
Laser Emission Indicator.....	20
Door Interlock.....	20
Emergency Stop	20
Protective Housing.....	20
Safety Interlocks	20
Safety Shutter	20
Location of Controls	20
Manual Reset.....	20
Electrical Fault Detection Circuitry	20
Location of Regulatory and Other System Labels	21
Ophthalmic Applications	26

Contraindications	26
Potential Complications or Side Effects	27
Adverse Effects and Complications	27
Posterior Segment Laser Procedures	27
Anterior Segment Laser Procedures	28
System Components	29
Laser Console	29
Touchscreen LCD Control Panel	29
Slit Lamp (model SL-PA05)	30
Micro Manipulator	30
Power Knob	30
Footswitch	31
Door Interlock Plug	31
Aperture Sheet	31
Breath Shield	31
Safety Glasses	31
Elbow Pad	32
Key	32
Power Cord	32
Laser Warning Sign	32
LIO (Laser Indirect Ophthalmoscope) – Optional Accessory	32
Wireless Footswitch – Optional Accessory	32
Endpoint Management Features (Software Option)	32
MicroPulse (Software Option)	33
References	34
Pattern Scanning Laser Trabeculoplasty (PSLT) (Software Option)	34
Indications for Use	34
Description	34
System Specifications	36
System Installation and Setup	39
Adjusting the Table Height	39
Connecting the System Components	40
Front Panel Connections	40
System Start-up and Shut-down	41
Starting the System	41
Shutting Down the System (Standard Shut-down)	41

Emergency Shut-down.....	41
Home Screen.....	42
System Setup Screen	43
Settings Tab.....	43
Patient Information Feature.....	44
Counter Reset.....	45
Treatment Time.....	45
Surgical Timeout	46
Physician Preferences	46
Advanced EPM Mode Display.....	46
Progressive Titrate	47
Audio Volume.....	47
Fire Volume.....	47
Treatment Report Window	47
Procedure to prepare the USB drive before the first time use	48
Posterior Reports	49
Anterior Reports	50
Auxiliary Control.....	50
Time Format.....	51
Date Format	51
Additional Features	51
Lens Tab	51
Posterior/Anterior Lenses.....	51
Custom Contact Lens Magnification	51
Physicians Tab.....	51
Auto-Transition Feature	53
Treatment Screen	53
Top Panel.....	53
Favorite Window	54
Edit Favorites Window	55
Center Panel	56
Controls.....	56
Power.....	56
Exposure.....	57
Spot Diameter Control	59
Spacing Field	60

Endpoint Management.....	62
MicroPulse	63
Radius Field	63
Curvature Field	65
Rep Rate Field	65
Treatment Plan Field.....	67
Titration	67
Internal Fixation Light.....	69
Outline.....	69
Auto Advance Feature	70
Auto-Transition Feature	71
AIM.....	71
Center	71
Rotation.....	72
Bottom Panel	72
Posterior Pattern Descriptions/Pattern Parameters	72
Single Pallet.....	73
Single Spot.....	73
Pattern Pallet	74
Array	74
Hexagon.....	74
Triple Arc.....	75
Single Spot.....	75
Wedge.....	75
Arc.....	76
Triple Ring.....	76
Line	77
Octant	78
Enhanced Octant	78
MicroPulse Pattern Pallet.....	79
EPM Pattern Pallet.....	80
Anterior Pallet	81
PSLT	81
Treatment Plan.....	81
Array	81
Single Spot.....	82

<i>Intraoperative Instructions</i>	83
Slit Lamp Treatment Procedure	83
LIO Treatment Procedure	84
Between Patient Treatments.....	84
System Shut-down	84
<i>Maintenance Instructions</i>	85
Annual Maintenance	85
System Repair	85
User Maintenance.....	85
Cleaning the Console External Surfaces	85
Cleaning the Control Panel Screen.....	85
Maintaining Effectiveness of Grounds	85
Replacing the Fuses	86
<i>Troubleshooting Guide</i>	86
System does not turn on	87
System monitor is blank for more than 30 seconds	87
Trouble Adding/Removing Favorites.....	87
No aiming beam is present when in READY mode and/or no laser treatment light is delivered when the footswitch is depressed and/or the beams are of poor quality	88
<i>Error Messages</i>	89
<i>Calibration Procedure</i>	100
Disclaimer Warning	100
Calibration Instructions	100
<i>System Relocation Instructions</i>	101
<i>Room Preparation</i>	101
<i>Instructions for Use</i>	101
Intended Use Environment	101
Patient Environment.....	102
<i>Electromagnetic Compatibility</i>	103
<i>Ophthalmology References</i>	108
<i>General information on Usage</i>	109
Intended Patient Population	109
Intended User Profile	109
<i>Warranty Information</i>	110
Warranty Shipments, Returns, and Adjustments	110

Disposal.....	111
Decontamination of Returned Equipment	112
Decontamination Certification	113
US Technical Service Information.....	113

Cautions and Warnings

Disclaimer

Calibration of the Iridex Pascal is a service procedure to be performed only by Iridex Corporation-certified personnel or customers who have taken and passed an Iridex Corporation Bio-Medical Preventative Maintenance Training course on the appropriate laser system. Adjustment by anyone other than Iridex Corporation-certified personnel or by customers who have taken and passed a Bio-Medical Preventative Maintenance Training course voids any existing manufacturer's warranty on the instrument and can result in serious personal injury.

Operator Manual Symbol Definitions

Please read this manual and follow its instructions carefully. The words "**WARNING**", "**CAUTION**", and "**NOTE**" carry special meanings and should be carefully reviewed.

	CAUTION	Alert the user to use special care necessary for the safe and effective use of the device. Cautions may include actions to be taken to avoid effects on patients or users that may not be potentially life threatening or result in serious injury, but about which the user should be aware. Cautions are also provided to alert the user to adverse effects on this device of use or misuse and the care necessary to avoid such effects.
	WARNING	Alert the user to potential serious outcomes (death, injury, or serious adverse events) to the patient or the user.
	LASER WARNING	Warning specifically related to hazard from a laser beam.
	DANGEROUS WARNING	Warning specifically related to a hazard from electricity.
	NOTE	Provided when additional general information is applicable.

General Safety and Regulatory Information

Iridex Corporation laser systems are precision medical instruments. The systems have undergone extensive testing. With proper handling, they are useful and reliable clinical instruments. To protect operating personnel and patients, this safety section and the appropriate slit lamp and pattern generator delivery system safety section should be read thoroughly before operation.

Iridex Corporation lasers are classified as Class IV lasers by the National Center for Devices and Radiological Health. Class IV represents the highest power lasers; for this reason, the user must take precautions to prevent exposure of laser energy to the eye and skin from either direct or diffusely reflected laser beams, except as a therapeutic application. In addition, precautions must be taken in the surgical environment to prevent the hazards of fire and electrical injury.

Iridex Corporation does not recommend specific clinical practices. The following precautions are extensive but may not be complete. Laser users are advised to supplement this information with technological advances in surgical products and techniques as they become available to the medical laser user community through medical literature. See also the American National Standards Institute (ANSI) publications.

- ANSI Z136.3-2005 – American National Standard for the Safe Use of Lasers in Health Care Facilities,
- ANSI Z136.1-2000 – American National Standard for the Safe Use of Lasers, CAN/CSA-S386-2008 – Laser Safety in Health Care Facilities and other national standards as may be applicable for the country in which the laser system is used.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: 1) This device may not cause harmful interference, and 2) this device must accept any interference received, including interference that may cause undesired operation.

Ocular Protection



WARNING: LASER HAZARD

Never look directly into the laser aperture or scattered laser light from reflective surfaces when the treatment beam is activated. Severe eye damage could occur.

Never look into the path of a laser beam. Laser safety eyewear only offers protection against stray or diffuse laser beam energy for a maximum exposure of 10 seconds.

Never substitute glass prescription eyewear for the appropriate laser safety eyewear, as severe eye damage could occur. The glass in prescription eyewear can concentrate the laser light onto the retina. A high power density beam can also shatter glass prescription eyewear, resulting in possible severe eye damage.

Do not use eyewear that is broken or damaged.

The  (laser emission) indicator is displayed on the Treatment screen to warn the user that the system is capable of emitting laser energy. Appropriate precautions, such as wearing appropriate eyewear in the room, should be taken.

As a precaution against accidental exposure to the output beam or its reflection, anyone checking or adjusting calibration should wear appropriate laser safety eyewear.

Because prolonged intense light exposure can damage the retina, the use of the device for ocular examination should not be unnecessarily prolonged and the brightness setting should not exceed what is needed to provide clear visualization of the target structures.

The retinal exposure dose for a photochemical hazard is a product of the radiance and the exposure time. If the value of radiance were reduced in half, twice the time would be needed to reach the maximum exposure limit.

While no acute optical radiation hazards have been identified for slit lamps, it is recommended that the intensity of light directed into the patient's eye be limited to the minimum level that is necessary for diagnosis. Infants, aphakes and persons with diseased eyes will be at greater risk. The risk may also increase if the person being examined has had any exposure with the same instrument or any other ophthalmic instrument using a visible light source during the previous 24 hours. This will apply particularly if the eye has been exposed to retinal photography. The red diode laser aiming beam has an average power varying from barely visible to 1 mW maximum. The safe (Class II) exposure duration limit at a maximum power level of 1 mW is 3.9 seconds. To protect the patient from possible retinal damage during treatment, use the lowest practical aiming beam intensity and the minimal required time duration.

Using any of the larger patterns while the slit lamp magnification is set to 32X, may result in the pattern overfilling the visual field. Reduce the slit lamp magnification or adjust the pattern size. Do not attempt treatment unless the entire pattern is visible.

Laser Safety Eyewear



WARNING: LASER HAZARD

Laser safety eyewear is routinely required with most lasers. When using the system, the laser safety officer should determine the need for safety eyewear based on the Maximum Permissible Exposure (MPE), Nominal Hazard Zone (NHZ) and the Nominal Ocular Hazard Distance (NOHD) for each of the available laser wavelengths, as well as the wavelength itself and the configuration of the treatment room (usually within the controlled area).

ANSI Standard Z136.1-2007 defines MPE as "the level of radiation to which a person may be exposed without hazardous effect or adverse biological changes in the eye or skin"; the NHZ as "the space within which the level of direct, reflected or scattered radiation during normal operation is not expected to exceed the applicable MPE"; and the NOHD as "the distance along the axis of the unobstructed beam from the laser to the human eye beyond which the irradiance or radiant exposure during operation is not expected to exceed the appropriate MPE."

The NOHD is measured from the slit lamp and pattern generator delivery system laser aperture. ANSI defines the controlled areas as "an area where the occupancy and activity of those within is subject to control and supervision for the purpose of protection from radiation hazards."

All personnel who are within the NOHD are considered to be in the controlled area and shall wear eye protection with the appropriate optical density. Eyewear must be resistant to physical damage and photo-bleaching. The minimum optical density (OD) is 4 at 532 nm or 577 nm; for countries inside Europe and that comply with EN 207, the eyewear must have a protection class of L5 at 532 nm or L4 at 577 nm.

Delivery Device	NOHD (532 nm and 577 nm)
Slit Lamp	5.4 m (17.7 feet)
LIO	16.1 m (52.8 feet)



NOTE

These ratings only apply to laser exposure greater than 200 mm (7.9 in) away from the laser exit aperture of the Slit Lamp.

The type of eye protection recommended for the physician, the patient, and/or treatment room personnel within the NHZ depends on the planned procedure and the equipment required to perform that procedure.

An eye safety filter is provided with the Slit Lamp and is required for safe use. Laser safety eyewear is not required by the physician who views the procedure through the slit lamp eyepieces. All other personnel within the NHZ must wear laser safety eyewear with the recommended optical density.

Along with providing the appropriate safety eyewear, the following steps should be taken to secure the controlled area:

1. Treatment should be conducted in a dedicated, enclosed room.
2. A warning sign should be placed on the outside of the treatment room door when the laser is in use. The sign is intended to alert personnel before they enter the controlled area.
3. The treatment room door should be kept closed during treatment.

Electrical Hazards



WARNING: DANGEROUS VOLTAGE

To avoid the risk of electric shock, this equipment must only be connected to a supply main with protective earth. Hospital Grade Cord grounding reliability can only be achieved when equipment is connected to an equivalent receptacle marked "Hospital Grade" or "Hospital only".

To avoid the risk of electric shock, do not touch any external connector and the patient simultaneously.

Do not use power cables other than the power cable provided with the system. Do not use extension cables with the system.

Disconnect the laser system from the electrical outlet when inspecting the fuses.

Never open the laser console protective covers. Opening the covers will expose you to high voltage components, the laser resonator and possible laser radiation. Only certified personnel shall work inside the console.

The area around the laser and footswitch shall be kept dry. Do not operate the laser if any of the cords are faulty or frayed. The laser should undergo routine inspection and maintenance per the Iridex Corporation manufacturer's recommendation and institutional standards.

Non-Homogeneous Pattern Delivery



WARNING

The optical system in this unit has been designed to provide uniform laser energy deposition over the entire area of a focused spot. If proper focus of the laser spot on its intended target is not achieved, or if damage, contamination, or deterioration of the optical components has occurred, it is also possible to degrade this uniformity. Observation of aiming laser spots should provide suitable indication of the level of uniformity attainable in any given system configuration, and treatment should proceed only when the user is satisfied as to the level of laser deposition uniformity by observation of properly focused aiming laser spots. Service personnel should be contacted for any concerns in this regard.

Treatment with laser deposition that is highly non-uniform can result in localized over- and/or under-treatment of the affected areas.

The larger the pattern, the more likely that spots within the pattern would be non-uniformly delivered. Smaller pattern sizes may be less likely to produce non-homogeneous lesions than larger patterns.

There are other variables which may contribute to non-uniform pattern laser application, including, but not limited to the following: media opacities (i.e. cataract) and heterogeneity within a particular opacity; ischemic retinal changes; other situations where there is no visible retinal/media non-homogeneity.

Fiber Optic Cable Assembly



WARNING: LASER HAZARD

Use extreme care with the cable assembly to/from the console. The cable assembly consists of wiring and fiber optic cables.

Do not pull or stress any cables. Do not exceed bend radius of 15 cm.

Do not set items on or under the cable assembly.

Damage to the fiber optic cables may cause unintended laser exposure.

Pattern Titration



WARNING

The physician is expected to show discretion in the dosage and location of the laser delivery when using the Pattern Titration feature. It is the responsibility of the physician to select the appropriate power and treatment location.

Please note that when retrieving a favorite that has been saved in titration mode, the values of the pattern parameters that are not allowed to be changed in titration will be restored to the defaults.

Titration



WARNING

A single titration session at the very beginning of a treatment may be insufficient, as different areas of tissue may respond differently over the course of the treatment. Uniform tissue response should be continuously assessed by the end-user and re-titration may be necessary during the procedure.

Please confirm that Endpoint Management is automatically disabled and grayed out during titration mode.

Attempting to perform titration when Endpoint Management is turned on may result in overpower exposure.

Fire Hazard



WARNING

Do not use the laser system in the presence of flammables or explosives such as volatile anesthetics, alcohol, certain surgical preparation solutions or other such substances. An explosion and/or fire could occur.

Do not use in an oxygen-rich environment.

The treatment beam can ignite most nonmetallic materials. Use fire retardant drapes and gowns. The area around the treatment site can be protected with towels or gauze sponges moistened with sterile saline solution or sterile water. If allowed to dry, protective towels and sponges can increase the potential fire hazard. A UL-approved fire extinguisher should be readily available.

Per IEC 60601-2-22: The use of flammable anesthetics or oxidizing gases such as nitrous oxide (N₂O) and oxygen should be avoided. Some materials (e.g., cotton wool) when saturated with oxygen may be ignited by the high temperatures produced in normal use of the system. The solvents of adhesives and flammable solutions used for cleaning and disinfecting should be allowed to evaporate before the laser system is used. Attention should also be drawn to the danger of ignition of endogenous gases.

Protecting Non-target Tissues



WARNING: LASER HAZARD

Never place hands or other objects in the path of the laser beam. Severe burns could occur.

Except during actual treatment, the system must always be in STANDBY mode. Maintaining the system in STANDBY mode prevents accidental laser exposure if the footswitch is inadvertently depressed.

Only the person aiming the laser beam should have access to the footswitch. Use caution depressing the footswitch when it is in proximity to a footswitch for other equipment. Make sure the footswitch depressed is the correct one to avoid accidental laser exposure.

Operational Safety



WARNING

Read this operating manual thoroughly and be familiar with its contents prior to using this equipment.

If excessive or unintended eye movement is a concern, treatment is not recommended.

Exercise caution while setting treatment parameters (for example exposure time and the number of spots per pattern) when laser burns are to be delivered in the non-macular area for long time periods leading to longer grid completion times. Please be aware that with longer completion times, the possibility of patient movement increases the risk of treatment of unintended targets.

Verify eye safety filter is properly attached to the slit lamp prior to use.

Always verify the power settings on the screen prior to pressing the footswitch.

Verify adjustments to laser parameters on the monitor before pressing the footswitch.

Verify that the slit lamp eyepieces are adjusted to your settings before each use, especially in a multi-user practice. Only when the eyepieces are properly adjusted is the laser confocal with the microscope. The laser spot diameter will not be accurate and may result in over treatment or under treatment if the eyepieces are not properly adjusted.

When the system is in READY mode, if the aiming beam is not present, is distorted, or is incomplete, do not proceed with treatment. Turn off the machine and contact service.

It is the responsibility of the physician to select appropriate combinations of repetition rate and exposure time to avoid overexposure or unintended exposure.

It is the responsibility of the physician to verify that the pattern visualized through the slit lamp is the same as the pattern displayed on the control panel. A discrepancy between patterns could indicate a hardware failure. Should this occur, discontinue treatment and contact service.

Early release of the footswitch will terminate the treatment beam before the complete pattern has been delivered. It is recommended that the pattern be completed by delivering the burns individually. DO NOT repeat the pattern at the same targeted tissue, as the pattern will start from the beginning, resulting in overlapping burns.

When a pattern with multiple spots is selected, use caution when operating with a multiple mirror contact lens. Do not overfill the mirror with the pattern and ensure that you have visualization of the complete pattern and the area to be treated prior to laser treatment.

Do not use any contact lens with a laser spot magnification of <0.94 .

Do not use wide field contact lenses. Wide field lenses will enlarge the spot diameter and alter the Foveal Exclusion Zone ring diameter.

Selecting the wrong contact lens or entering the incorrect magnification factor for a custom lens will cause the displayed fluence to be incorrect. Endpoint Management modulates both Power (mW) and Exposure (ms) simultaneously. When Power reaches the lower limit, only Exposure is adjusted. If engaging Endpoint Management at or near these minimum power limits, lower Endpoint Management percentage values may be required.

Do not use the LIO if it has been dropped or if there is visible damage to the housing or fiber until the LIO has been checked to work properly by certified personnel.

The video monitor must not be used to guide treatment or for diagnosis. The treating physician must visualize treatment through the slit lamp at all times.

Footswitch marked with IPx8 is suitable for general or office use. Do not use the system with footswitch in the operating room.

When used in a surgical environment, ensure all concentrations of O₂ are low and that flammables or volatile anesthetics, alcohol and surgical preparation solutions are decreased. All solution-soaked materials must be separately bagged and contained and/or removed from the room prior to the use of laser (e.g. prepping dispensers, applicators, drip drapes). This includes items used by the anesthesia providers. The anesthesia provider will decrease FIO₂ (Fraction of inspired oxygen) to minimum amount possible to maintain adequate SpO₂. Be aware of possible enriched O₂ and N₂O atmospheres near the surgical site under the drapes, especially during head and neck surgery. Tent drapes to allow oxygen, which is slightly heavier than air to drain away from the patient's head and toward the floor.



CAUTION

If the system becomes unresponsive at any time other than during laser emission, do not press the emergency laser stop button. Instead, turn the key switch to the OFF position. Wait at least one minute before restarting the system using the key switch.

If the control panel is blank for more than 60 seconds during system start-up, verify that the power indicator LED on the front of the control panel is illuminated. If it is not illuminated, then press the Power button on the right of the control panel to turn on the control panel. If the control panel remains blank, turn off the system with the key; verify that all control panel cables are plugged in and fully seated; and then restart the system. If the screen is still blank, turn off the system and contact service.

It may take longer for the equipment to achieve a ready state in a low-temperature environment.

It is the responsibility of the physician to verify that the aiming beam spot visualized through the slit lamp is the expected size. If the aiming beam size or pattern appears inappropriate or distorted, do not proceed with treatment. Readjust the slit lamp focus. If the problem persists, contact service.

It is the responsibility of the physician to select the appropriate power and treatment location. The lowest practical setting should always be used to achieve the desired clinical outcome.

Do not use a wet cloth to clean the control panel screen. Doing so may damage the screen.

Additional Safety Considerations



WARNING

US federal law restricts this device to sale by or on order of a physician. (CFR 801.109(b)(1))

Use of controls or adjustments or performance of procedures other than those specified herein may result in hazardous radiation exposure.

Iridex Corporation medical devices are solely for use by physicians trained in the operation of laser photocoagulation and associated delivery devices.

To avoid potential injury to the user and the patient and/or damage to this device, the user must:

- *Read this manual thoroughly and be familiar with its contents prior to using this equipment.*
- *Be a qualified physician, having complete knowledge of the use of this device.*
- *Test this device prior to a procedure.*
- *Attempt no internal repairs or adjustments not specifically detailed in this manual.*

Do not modify this equipment without authorization of the manufacturer.

When the laser system is interconnected with other medical electrical equipment, leakage currents may be additive. Ensure all systems are installed according to the requirements of IEC 60601-1.

If the laser system is used adjacent to or stacked with other equipment, observe and verify normal operation of the laser system in the configuration in which it will be used prior to use.

Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the laser system, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.



CAUTION

Condensation may occur if the equipment is exposed to high humidity for an extended period of time.

Vibration or physical shock may affect the quality, performance and reliability of the equipment.

Do not connect two footswitches to the laser console.

Safety Compliance

Caution: US Federal law restricts this device to sale by or on the order of a physician.

The Iridex PASCAL® (532 nm and 577 nm models) conform to the clauses of IEC 60825-1 Ed. 3 and IEC 60601-2-22 Ed. 4 identified in Section III of Laser Products – Conformance with IEC 60825-1 Ed. 3 and IEC 60601-2-22 Ed. 4 (Laser Notice No. 56) as comparable with 21 CFR 1040.10 and 1040.11.

The UDI Product Iridex PASCAL Label, provided on the device, includes:

- Identification information as specified in 21 CFR 1010.3
- The statement “Complies with FDA performance standards for laser products except for conformance with IEC 60825-1 Ed. 3 and IEC 60601-2-22 Ed. 4, as described in Laser Notice No. 56, dated May 8, 2019.” to meet the certification requirements of 21 CFR 1010.2
- The Prescription Only symbol, as specified in 21 CFR 801.109

Regulatory Compliance Safety Features

Iridex Pascal complies with 21 CFR subchapter J as administrated by the Center for Devices and Radiological Health of the Food and Drug Administration (FDA). The following FDA compliance safety features are included:

Key Lock Switch

The system can be activated only with the proper key to operate the master key switch. The key cannot be removed in the ON position and the system will operate only with the key in place. When treatments are complete, always remove and secure the key to prevent unauthorized use of the system.

Laser Emission Indicator

The laser emission indicator is displayed to warn the user that the system is capable of emitting laser energy and that appropriate precautions should be taken, such as using the appropriate eyewear when in the treatment room.

Door Interlock

A door interlock may be used in conjunction with a remote switch to disable the system in case of certain external events (e.g., the opening of a treatment room door). A remote switch or interlock can be wired to the door interlock plug and connected to the system interlock receptacle on the rear of the system console. If a remote switch is used, the system can be set in READY mode only when the remote switch is closed. Breaking the connection by opening the switch (door) or removing the plug disables the system and the system returns to STANDBY mode with "Door Interlock" displayed on the control panel.

Emergency Stop

When pressed, immediately turns off power to the laser.

Protective Housing

The system console has a protective housing that prevents unintended human access to laser radiation above Class I limits. This housing is to be opened only by certified personnel.

Safety Interlocks

The protective housing is not designed to be removed by the user during operation or maintenance. Therefore, the system does not have and is not required to have any safety interlock within the meaning of US FDA 21 CFR, Section 1040 or European EN 60825-1.

Safety Shutter

The laser system utilizes an electronic laser safety shutter. The system will not be able to emit laser light unless all safety conditions are met prior to depressing the footswitch. The safety shutter is activated when the system is off, during the self-test at power-on, in STANDBY mode, or when the safety monitor detects a fault.

Location of Controls

Controls are located on the touchscreen control panel.

Manual Reset

If laser emission is externally interrupted during treatment by activation of the door interlock, the system will automatically go into STANDBY and the safety shutter will revert to a closed position. To resume treatment, reset the system by placing the laser in READY mode.

If laser emission is interrupted by main electrical power loss, the system will automatically turn off. To resume treatment after an electrical power loss, the system must first be manually restarted by rotating the key switch to the ON position.

Electrical Fault Detection Circuitry

If the electronic system detects a fault condition, laser exposure cannot occur. The laser is disabled,


Laser Radiation Caution Label, includes:



Iridex Pascal (577nm)




Laser Emission Warning

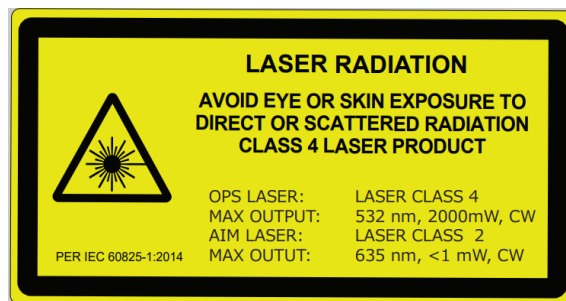

Iridex Pascal (532 nm)

Caution, Class 4 Laser Emission


Laser beam may be present in this area, when opened.

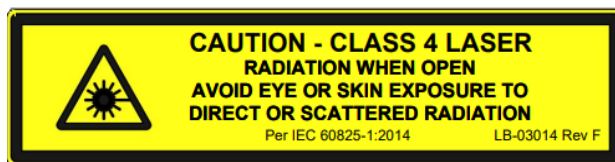

Laser Aperture Label

Laser Emission Warning




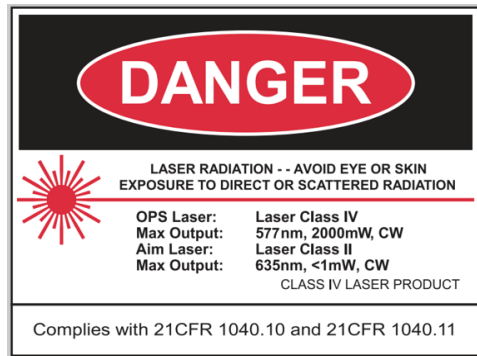


Iridex Pascal (577nm)

Danger Label, includes:

- Laser Emission Warning
- Wavelength
- Power
- Laser Class

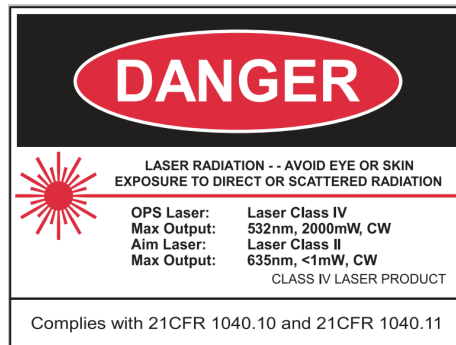

Iridex Pascal (532 nm)

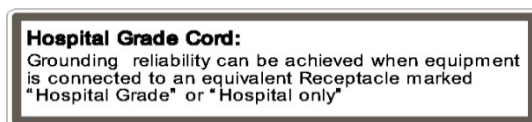

Hospital Grade Cord Grounding Instructions



USB Connections



Remote Door Interlock Connection



Footswitch Connection


I

ON

O

OFF

USB-C Cable for Monitor

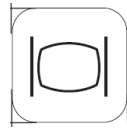

LIO Connection



IRIDEX® Iridex PASCAL® 532

Iridex Corporation
1212 Terra Bella Avenue
Mountain View, CA 94043 USA
+1.650.940.4700
pascalservice@iridex.com
iridexretina.com

100 - 240 ~, 50/60 Hz
600 VA
Tested and Complies with FCC Part 15 Class B

Complies with FDA performance standards for laser products except for conformance with IEC 60825-1 Ed.3 and IEC 60601-2-22 Ed. 4, as described in Laser Notice No. 56, dated May 8, 2019

Conforms to:
IEC 60601-1:2005+A1:2012+A2:2020
IEC 60825-1:2014
IEC 60601-2-22:2019
IEC 60601-1-6:2010 + A1:2013 + A2:2020
ANSI/AAMI ES60601-1:2005 +A2: 2021
CAN/CSA-C22.2 No. 60601-1:14

REF SA-06567

SN 123456

W 2026-02-18



(01)00813125016890
(11)260218
(21)123456

The instructions for use and symbols glossary are provided in the Operator Manuals



<170lbs. (<70kg)

R_X ONLY



LB-06131 Rev E

System Information Label, includes:



Manufacturer



Part Number



Serial Number



Manufacture Date



Follow operating instructions



WEEE Directive



Type B applied part



Caution



Prescription Only



Medical Device



Unique Device Identifier



Space to apply the following ETL mark:



ETL mark



Mass Symbol



Fuse

IRIDEX® Iridex PASCAL® 577

Iridex Corporation
1212 Terra Bella Avenue
Mountain View, CA 94043 USA
+1.650.940.4700
pascalservice@iridex.com
iridexretina.com

100 - 240 ~, 50/60 Hz
600 VA
Tested and Complies with FCC Part 15 Class B

Complies with FDA performance standards for laser products except for conformance with IEC 60825-1 Ed.3 and IEC 60601-2-22 Ed. 4, as described in Laser Notice No. 56, dated May 8, 2019

Conforms to:
IEC 60601-1:2005+A1:2012+A2:2020
IEC 60825-1:2014
IEC 60601-2-22:2019
IEC 60601-1-6:2010 + A1:2013 + A2:2020
ANSI/AAMI ES60601-1:2005 +A2: 2021
CAN/CSA-C22.2 No. 60601-1:14

REF SA-06568

SN 123456

W 2026-02-18



(01)00813125016906
(11)260218
(21)123456

The instructions for use and symbols glossary are provided in the Operator Manuals



<170lbs. (<70kg)

R_X ONLY



LB-06131 Rev E



Footswitch	Ingress Protection
IPX 8	Suitable for continuous immersion water
	Caution
CE 2797	CE mark with BSI Netherlands NB Number

Ophthalmic Applications

Iridex Pascal is used to treat ocular diseases in both the posterior and anterior chambers of the eye. The laser systems are well-suited for treating the eye because they have minimal effect on transparent tissues and materials. Iridex Pascal laser energy can be efficiently delivered to opaque structures of the eye through the transparent cornea, aqueous humor, lens and vitreous humor, allowing many conditions to be treated by non-invasive techniques. Laser energy is delivered to opaque structures within the eye through a slit lamp that has been specially adapted for use as a laser delivery system.

The delivery system includes a lens system to focus the laser energy and vary the size of the laser spot in the plane of observation of the slit lamp. It includes a mechanism to manipulate the position of the laser beam without moving the slit lamp. Laser energy is delivered to the slit lamp through a flexible fiber optic.

For most procedures, a laser contact lens is used to direct laser energy to the part of the eye being treated. The contact lens may have mirrors so that laser energy can be delivered to areas of the retina behind the iris or into the angle so that the trabecular meshwork can be treated. The contact lens also helps to hold the eye open and still so that laser energy can be delivered effectively.

Systems may be used for procedures performed in a hospital or in a physician's office for in-patient or out-patient procedures. The use of the laser system is not a contributing factor in deciding whether a procedure is done on an in- or out-patient basis.

Contraindications

The following conditions are contraindications for performing laser trabeculoplasty:

- Any corneal opacities, cataract formation and vitreous hemorrhage which can interfere with the laser surgeon's view of the target structures
- Aphakic eye with vitreous in the anterior chamber
- Neovascular glaucoma
- Glaucoma caused by congenital abnormalities of the angle
- Less than 90° of open angle or extensive low-lying peripheral anterior synechiae present circumferentially around the angle
- Significant corneal edema or a diminished aqueous clarity obscuring visualization of the angle detail
- Glaucoma secondary to active uveitis

Tissue absorption is directly dependent upon the presence of pigmentation; therefore, dark pigmented eyes will require lower energies to obtain equivalent results as compared to light pigmented eyes. Do not treat albino patients who have no pigmentation.

For patients with wide variations in retinal pigmentation as evaluated by ophthalmoscopic observation, select multi-spot patterns which cover a homogeneously pigmented smaller area to avoid

unpredictable tissue damage.

Exercise caution while setting treatment parameters (for example exposure time and the number of spots per pattern) when laser burns are to be delivered in the non-macular area for long time periods leading to longer grid completion times. Please be aware that with longer completion times, the possibility of patient movement increases the risk of treatment of unintended targets.

In addition, following are the contraindications for photocoagulation treatment:

- Patient cannot fixate their eye or hold still (for example patient has nystagmus).
- There is inadequate view of the fundus due to opacity (not clear media so physician cannot see the fundus).
- Presence of subretinal fluid in the patient's eye.

Potential Complications or Side Effects

Potential complications specific to retinal photocoagulation include inadvertent foveal burns, choroidal neovascularization, paracentral scotomata, subretinal fibrosis, photocoagulation scar expansion, Bruch's membrane rupture, choroidal detachment, exudative retinal detachment, pupillary abnormalities from damage to the ciliary nerves, and optic neuritis from treatment directly or adjacent to the disc.

Potential complications specific to laser iridotomy include iritis, visual symptoms, and rarely, retinal detachment.

Adverse Effects and Complications

Posterior Segment Laser Procedures



WARNING

The most common complication of panretinal photocoagulation is increased macular edema usually with a concurrent decrease in visual acuity. In addition, blowout hemorrhages from the areas of neovascularization, particularly on the optic nerve, have been observed and may be caused by an increase in peripheral resistance secondary to photocoagulation or by an inadvertent valsalva maneuver by the patient.

Only a contact lens specifically designed for use with laser energy should be used. Use of a standard diagnostic contact lens may result in a power loss due to reflection from the surface of the lens. The reflected energy may pose a hazard to both the patient and the physician.



WARNING

Following photocoagulation, patients should be cautioned against any activity that could increase the venous pressure in the head, neck or eyes, such as straining, lifting or holding their breath. Patients should be advised to sleep with the head of their bed elevated 15 to 20 degrees.

Patients should be cautioned against stifling a sneeze, because this raises the blood pressure within the eyes to a high level. Vigorous nose blowing should also be discouraged. Rubbing the eyes following photocoagulation may disrupt blood vessels inside the eyes. Sneezing and coughing should be controlled with cough syrup or other medications.

Immediately following treatment, patients should avoid altitudes over 2500 m (~8000 ft.).

Anterior Segment Laser Procedures



WARNING

Intraocular pressure should be closely monitored following laser iridotomy or trabeculoplasty.

Hemorrhage from the trabecular meshwork occasionally occurs as an ooze of blood from Schlemm's canal to the site of laser impact. This is easily stopped by increasing the pressure on the gonio lens on the cornea or by coagulating the bleeding site by application of a laser burn.

Pupillary distortion may be encountered if the iris root or peripheral iris has been treated. This distortion may or may not be permanent, depending on the severity of the accidental damage.



WARNING

Intraocular pressure elevations have been reported to occur in up to 53% of eyes when 360° of the trabecular meshwork has been treated with 100 spots at the initial session. Intraocular pressure rises occur most frequently from 1 to 2 hours following laser treatment, although they may occur several hours afterward. For this reason, it is imperative to monitor patient intraocular pressure after laser treatment for up to 24 hours.

Peripheral anterior synechiae may occur when the posterior portion of the trabecular meshwork or other structures posterior to the meshwork are treated. These are best avoided by meticulous delivery of a well-focused laser beam.

Transient corneal epithelial burns have reportedly been resolved within 1 week without scarring. Endothelial burns are rarely encountered when careful focusing is employed.

Rarely, severe iritis may occur, related to either an unusual patient response or improper spot location.

System Components

Laser Console


Laser system under tabletop houses the key switch, emergency laser stop button, fiber ports, laser fiber, and electronics cable, 532 nm or 577 nm single wavelength laser, control electronics and power supply.

Laser Console




Touchscreen LCD Control Panel

Provides controls for selecting treatment parameters and displays for monitoring system information.

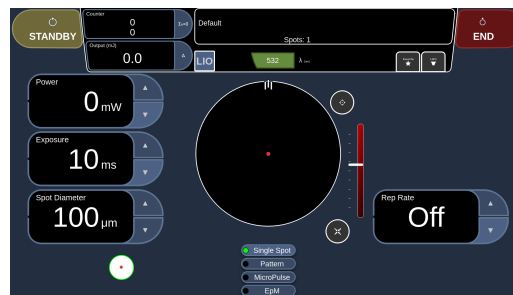

Slit Lamp (model SL-PA05)

Integrated slit lamp connection. Service personnel will connect the laser system to the slit lamp.

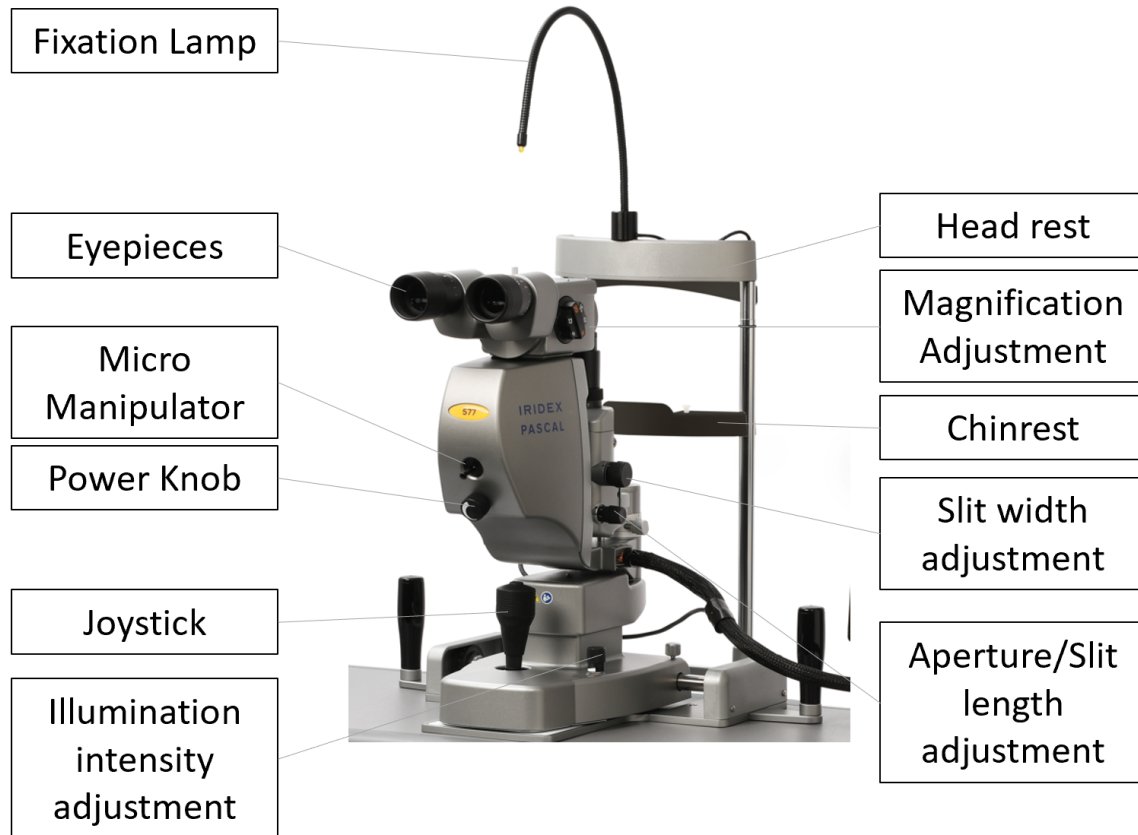

Iridex Pascal Slit Lamp PA05

Micro Manipulator

Provides auxiliary input for laser position. The aiming beam/treatment beam can be moved up/down/left/right in same direction as the Micro Manipulator (MM) moves.

The Micro Manipulator will return to mechanical center when released. There may be a few counts from true center after releasing the Micro Manipulator. Select the Home button to align the Micro Manipulator to center.



NOTE

The Micro Manipulator is disabled during the laser treatment.

Power Knob

The Power Knob located below the Micro Manipulator provides auxiliary input for laser power.

You can adjust power up or down by rotating the power knob similar to the touchscreen power up or down. Clockwise rotation increases the power and counter-clockwise rotation decreases the power.



NOTE

The Power Knob is disabled during the laser treatment.

Footswitch

Activates laser treatment beam when depressed while the system is in READY mode.



NOTE

Footswitch marked with IPX8 is suitable for general or office use.


Door Interlock Plug

Disables laser if the treatment room doors are opened or if the interlock plug is removed. Use is optional; however, the interlock plug must be inserted for the laser to operate.

Aperture Sheet

Provides an option to reduce the amount of reflections from the slit lamp.

Insert an aperture sheet as shown below.

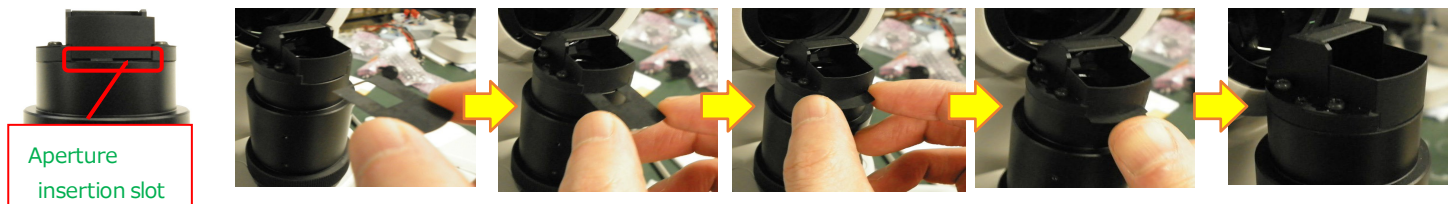

NOTE


Be sure to hold the arc-shaped side when inserting.

You can insert the sheet either side up.

Make sure to insert the sheet all the way as shown in the last picture above.

After inserting the sheet, close with the sticker provided (see photo of sticker on the right).

Exercise caution when handling the aperture sheet. The edges are sharp!

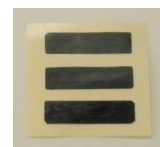

Breath Shield

Fits onto the slit lamp and provides a first line of defense for both healthcare providers and patients.

Safety Glasses

Provides eye protection while worn in the treatment room.

Elbow Pad

Stackable pad to use as an elbow rest.

Key

Power Cord

Laser Warning Sign

LIO (Laser Indirect Ophthalmoscope) – Optional Accessory

The LIO is an optional accessory to the Iridex Pascal system. Refer to the LIO Operator Manual for detailed safety and regulatory compliance information.

Wireless Footswitch – Optional Accessory

Endpoint Management Features (Software Option)

EpM (Endpoint Management): Refers to a software feature that allows delivery of laser energy which is a user-selected percentage of a laser dose pre-determined by titration. This pre-determined titration dose is the “100%” level, corresponding to the energy produced by the displayed “power” and “exposure” parameters on the Treatment screen, and is expected to be determined by ophthalmoscopically visible (sub-visible) test burns delivered by the physician. With Endpoint Management enabled, the laser power and exposure duration are reduced to deliver the user-selected percentage of the output energy settings. For example, if the user titrates to a power/exposure duration setting that delivers 4 mJ laser energy, an EpM setting of 75% would provide 3 mJ for each exposure. EpM percentage ranges from 10% to 95%.

Duty cycle within each laser pulse remains constant (100% duty cycle) with EpM.

To increase the dose above the 100% level, the user is expected to re-titrate laser power with additional test burns. By providing fine control over delivered energy, EpM allows the user to control the laser output to levels where ophthalmoscopically visible lesions are not achieved, while referencing a dosage with a visible effect (the 100% dosage).

The EpM features are only enabled on retinal treatment patterns for use in retinal photocoagulation. The use of EpM is otherwise restricted by the physician.

LM (Landmark): Refers to a feature of the Endpoint Management software application that allows delivery of two energy dose levels within a single pattern. When EpM is enabled, the user can optionally enable Landmark patterns. With EpM enabled, the outermost spots, the LM exposures, in the pattern are set to the 100% dosage (100% nominal power and exposure duration displayed on the Treatment screen), while the inner spots are delivered at the current EpM % setting. With EpM enabled and Landmarks disabled, the entire pattern is delivered at the current EpM % setting.


The purpose of LM patterns is twofold – to indicate the location of patterns delivered with EpM, which may be less ophthalmoscopically visible than the 100% exposures, and to provide visual feedback to the physician for dosimetry. In treatments without the use of Endpoint Management, it is routine for physicians to use the visual appearance of lesions to guide their adjustment of laser power in order to maintain a constant lesion grade.

Exposures with a low EpM % setting do not provide such guidance, but by delivering the outermost

spots (the Landmarks) in the pattern at the full 100% titration dose, this visual cue is maintained. As the Landmark burns in the pattern vary in effect as the laser application is moved across the retina, the user can adjust the laser power to maintain the same lesion appearance as the original titration burn.

MicroPulse (Software Option)

MicroPulse (μ P) is a laser delivery consisting of a group of microsecond bursts.

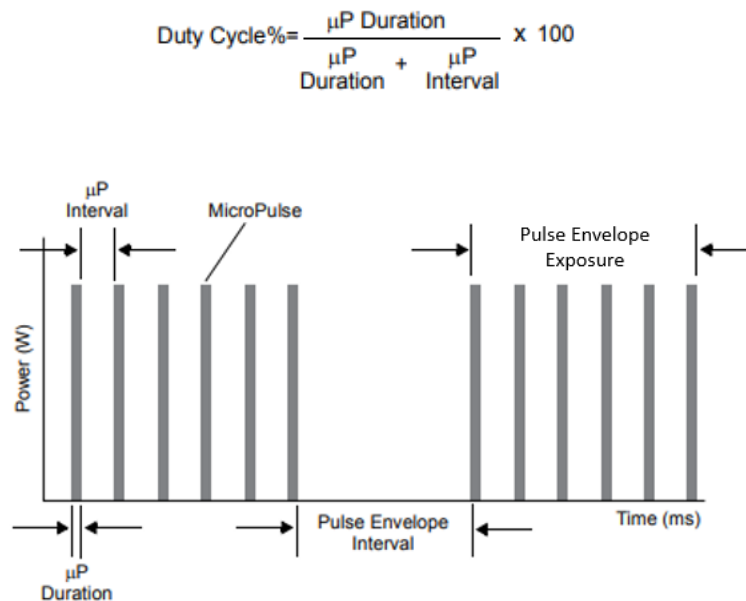

MicroPulse laser is typically used to administer subvisible threshold laser treatments to macular and perimacular targets. When used here, the terms “subvisible”, “subvisible threshold” or “subthreshold” denote that the desired endpoint is one in which treated tissue offers no ophthalmoscopically observable laser effects. Nevertheless, 577 nm and 810 nm studies have confirmed that subvisible laser treatment strategies can be clinically effective while inducing no changes discernible by slit lamp observation, fluorescein angiography (FA), fundus autofluorescence (FAF), or at any time postoperatively.^{1,2}

Tissues receiving subvisible MicroPulse laser treatment show no such changes because:

- MicroPulse laser delivery is being used instead of CW, and
- The total laser energy of such doses is only a percentage (often chosen by clinicians to be 20-70%) of that energy needed to produce a visible endpoint.

Energy (J) is equal to [Laser Power (W)] × [Exposure Duration(s)] × [Duty Factor (%/100)]. Duty Factor is often 5% to 15% when using MicroPulse mode and is 100% when using CW mode. Clinicians have reported various strategies to adjust these parameters relative to suprathreshold burns in order to achieve clinically effective subvisible endpoints.¹⁻⁴

Additional parameters to consider in any laser treatment protocol, and particularly during MicroPulse, is spacing between laser treatment spots, and the total number of treatment spots administered. Due to the limited thermal spread of MicroPulse exposures, subvisible treatments often call for the administration of a greater number of treatment spots with denser spacing than that used for threshold laser grid treatments.⁴

References

- ¹ Vujosevic S, Bottega E, Casciano M, Pilotto E, Convento E, Midena E. Microperimetry and fundus autofluorescence in diabetic macular edema: Subthreshold micropulse diode laser versus modified early treatment diabetic retinopathy study laser photocoagulation. *Retina* 2010;30(6):908-916.
- ² Vujosevic S, Martini F, Convento E, Longhin E, Kotsafti O, Parrozzani R, Midena E: Subthreshold Laser Therapy for Diabetic Macular Edema: Metabolic and Safety Issues. *Curr Med Chem* 2013;20(26):3267-71.
- ³ Figueira J, Khan J, Nunes S, Sivaprasad S, Rosa A, de Abreu JF, Cunha-Vaz JG, Chong NV. Prospective randomised controlled trial comparing sub-threshold micropulse diode laser photocoagulation and conventional green laser for clinically significant diabetic macular oedema. *Br J Ophthalmol* 2009;93(10):1341-4.
- ⁴ Lavinsky D, Cardillo JA, Melo LA, Jr., Dare A, Farah ME, Belfort R, Jr. Randomized clinical trial evaluating mETDRS versus normal or high-density micropulse photocoagulation for diabetic macular edema. *Invest Ophthalmol Vis Sci* 52(7):4314-23.

Application	Spot Size at Target* (µm)	Power (mW)	Duty Cycle (500 Hz)	Exposure Duration (ms)
Central Retina Focal/Grid	50–200	100–400	5%, 10%, 15%	100–300
Trabeculoplasty	200–300	400–1200	5%, 10%, 15%	100–300
*Spot size at target is dependent on many parameters, including spot size selection, physician's choice of laser delivery lens, and patient's refractive power.				

Pattern Scanning Laser Trabeculoplasty (PSLT) (Software Option)

Indications for Use

Pattern Scanning Laser Trabeculoplasty (PSLT) software is intended for use with the Iridex Pascal for procedures in trabeculoplasty in Open Angle Glaucoma.

Description

Pattern Scanning Laser Trabeculoplasty or PSLT is an advanced tissue-sparing laser treatment for reducing intraocular pressure. PSLT provides a rapid, precise, and minimally traumatic (sub-visible) computer-guided treatment that applies a sequence of patterns onto the trabecular meshwork. Automated rotation of consecutive patterns ensures that treatment steps are precisely placed around the trabecular meshwork without overlap or excessive gaps.

When the PSLT pattern is selected, the pattern is delivered when you depress the footswitch, after which the system automatically rotates the pattern clockwise. If the 360-degree treatment plan is selected, the pattern rotates a full 360 degrees before treatment is automatically terminated. If the 180-degree treatment plan is selected, the pattern rotates 180 degrees before treatment is automatically terminated.

1. Establish starting location for PSLT pattern.

2. Titrate to set desired Power.
3. Align to trabecular meshwork using Curvature.
4. "Treatment Complete" displays in the PSLT Treatment window, and the system enters STANDBY mode. Select **OK** to return to the Anterior Treatment screen.

Additional treatments can be initiated after returning to the Anterior Treatment screen.

Exposure Time	Spot Diameter	Curvature	Treatment Plan	Power
5 ms (Titrate off); 10 ms (Titrate on)	100 µm	0.00-3.00	360° or 180°	0-1500 mW



NOTE

For use with a gonioscopic lens of 1.0X magnification.

The PSLT pattern appears in the center of the field of view and not in the periphery. Use contact lens to set location.

If, during the course of the treatment, the footswitch is released before the entire pattern is delivered, "Treatment Complete" and "<NN Incomplete Segments>" (where "<NN>" represents the number of segments that were not completed) display in the PSLT Treatment summary window.

System Specifications

[Specifications are subject to change without notice.]

Treatment Beam		
	Iridex Pascal (532 nm)	Iridex Pascal (577 nm)
Type	Optically Pumped Semiconductor Laser (OPSL)	
Wavelength (nm)	532	577
Power output (mW)	0 – 2000 Power increments • 10 mW (from 30 mW to 150 mW, omitting 10 mW and 20 mW) • 25 mW (from 150 mW to 2000 mW)	
Duty cycle	100%	
Pulse durations (ms)	• Standard Package: 10 to 1000ms • PSLT (Optional): 5 to 1000ms • MicroPulse (Optional): 0.05 ms – 1.0 ms	
Pulse interval	1, 1.5, 2, 3, 4, 5, 6, 7 and 8 Hz (single spot or LIO)	
Pulse counter	0 – 99,999	
Laser beam diameter	50, 100, 200, 300 µm (in air)	
CDRH classification	Class IV	
European MDD laser classification	Class 4	
Aim Beam		
Type	Laser Diode	
Wavelength (nm)	635	
Power output	< 1 mW	
CDRH classification	Class II	
European MDD laser classification	Class 2	

Electrical Requirements		
	Iridex Pascal (532 nm)	Iridex Pascal (577 nm)
Voltage	100-240V~, 50/60Hz	
Rated power	600VA	
Fusing	218 Series, 5×20 mm, Time-Lag Fuse	
Fan noise	< 55 dBA	
Product Classifications per IEC 60601-1		
Class I equipment		
Type B equipment		
Standard equipment, Footswitch is IPX8		
Non-sterile product		
Equipment not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide		
Continuous operation		
Classifications & Approvals		
EN/IEC 60601-2-22	Laser Safety Requirements for Diagnostic and Therapeutic Laser Equipment	
EN/IEC 60601-1	International Safety Requirements for Medical Electrical Equipment	
EN/IEC 60601-1-2	EMC Requirements for Medical Electrical Equipment	
ISO 14971	Risk Management for Medical Devices	
CAN/CSA-C22.2 No. 60601-1	Canadian deviations for Medical Electrical Equipment	
ANSI/AAMI ES 60601-1	US Safety Requirements for Medical Electrical Equipment	
JIS T 0601-1	Japanese Requirements for Medical Electrical Equipment	
EN/IEC 60825-1	Safety of laser products	
FCC	Tested and Complies With FCC Part 15 Class B	

Environmental Requirements (operating)		
	Iridex Pascal (532 nm)	Iridex Pascal (577 nm)
Maximum altitude	2,000m (6562 ft.)	
Operating temperature	15° – 35°C (59° – 95°F)	
Maximum humidity	15° to 25°C: 85% (non-condensing) 25° to 35°C: 60% (non-condensing)	
Atmospheric pressure range	80.0 – 106.0 kPa	
Environmental Requirements (non-operating)		
Maximum altitude	Standard commercial shipping altitude	
Non-operating temperature	-10° to +55°C (14° to 131°F)	
Maximum humidity	85% (non-condensing)	
Atmospheric pressure range	70.0 – 106.0 kPa	
Physical Characteristics		
	Iridex Pascal (532 nm)	Iridex Pascal (577 nm)
Table height with casters	29 in (74 cm)	
Table width	36 in (92 cm)	
Table depth	18 in (46 cm)	
System weight	< 170 lbs. (< 77 kg)	
Power cable length	15 ft. (4.6 m)	
Laser fiber and electronics cable	Fiber bundle: 2.7m	
Footswitch cable length	4.5 m (180 in)	
Latex	This product is latex free	
Laser Safety Eyewear		
	Iridex Pascal (532 nm)	Iridex Pascal (577 nm)
Non-CE eyewear	Minimum OD of 3.8 at 532 nm per ANSI Z136.1	Minimum OD of 3.8 at 577 nm per ANSI Z136.1
CE eyewear	L5 at 532 nm per EN 207 Personal Eye Protection	L4 (minimum OD of 4) at 577 nm per EN 207 Personal Eye Protection

System Installation and Setup

Iridex Pascal is designed for installation and use in a darkened office or surgery room. Optimal system performance and viewing are achieved in low ambient light conditions. Iridex Pascal system installation and testing will be performed at your facility by Iridex Corporation-certified personnel.

Iridex Pascal is equipped with a Hospital Grade AC power cord. When selecting the location for system installation, ensure that the AC wall power outlet is correctly grounded. Follow local electrical codes to ensure proper grounding of the AC wall power outlet. A correctly grounded power connection is required for safe system operation.



NOTE

Do not position the laser system to make it difficult to operate the disconnection device.

Select an appropriate location that can accommodate the system size and allow easy access by both the patient and the physician. Ensure proper ventilation, temperature, and relative humidity. Select a well-ventilated space in an office or surgery room. See the “**Environmental Requirements (operating)**” section and ensure that the location of installation meets the temperature and relative humidity requirements listed. Position the system to aim the treatment beam away from windows and doors. Post a laser safety sign at the entrance to the treatment room.

Do not block cooling airflow or cooling vents on the laser system. Allow at least 5 cm (2 in) of clear space around the laser system to provide adequate system cooling airflow. Use care when routing system cables to prevent a tripping hazard and to protect the optical fiber from damage by being crushed under foot or being rolled over by a chair. If the cord must cross a floor where there is traffic, use of a floor cord/cable cover is recommended.


Adjusting the Table Height

Table height can be adjusted using the up/down arrow control as shown below.




Connecting the System Components

Refer to the following diagrams for the location of system component connections on the front and rear panels of the console.

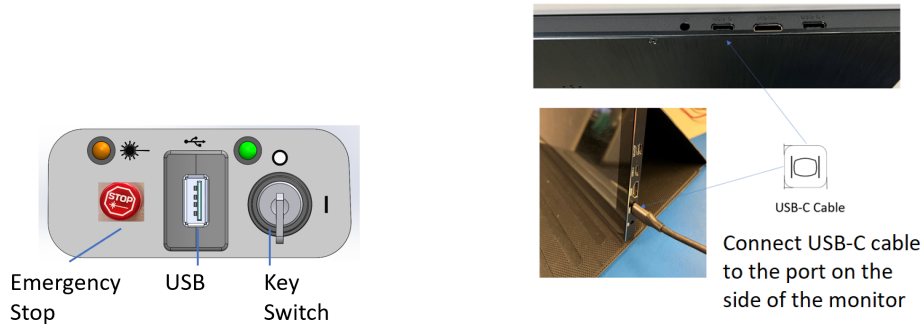

Front Panel Connections

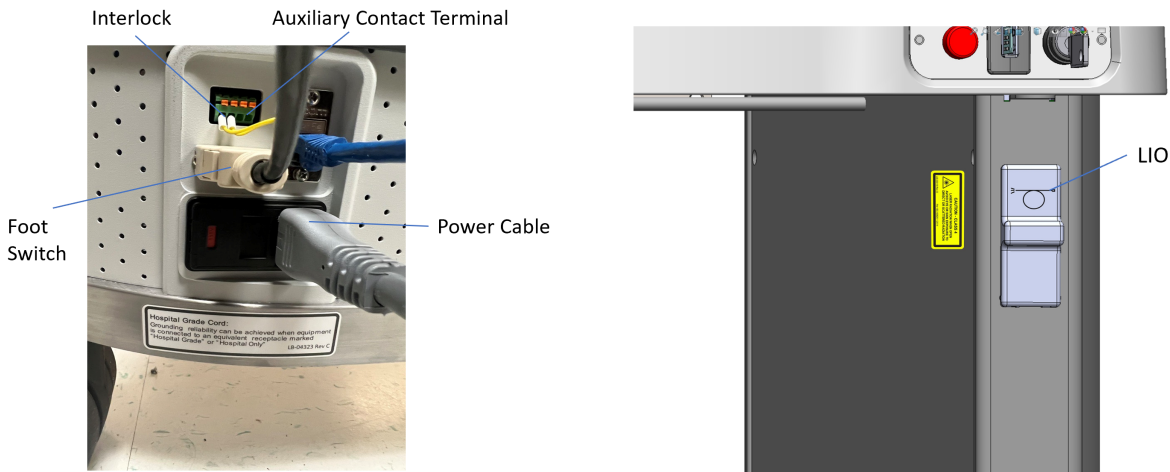

NOTE

If using an external door interlock, a qualified electrical professional must install the external switch, and the total length of the cable should not exceed 5 m (16 ft.).




CAUTION

USB Ports are not compatible with wireless devices. Use of the USB Ports with a wireless device cannot guarantee the performance of the wireless device or the Iridex Pascal System.

System Start-up and Shut-down

Starting the System

1. Connect the system to a wall power outlet.
2. Insert the key into the key switch.
3. Turn the key switch to the ON position.

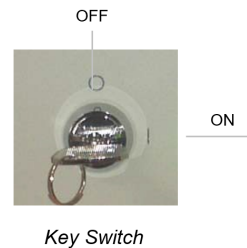

NOTE

Footswitch must be plugged in before powering up otherwise error will pop up.

Shutting Down the System (Standard Shut-down)

From the Treatment screen:

1. Place system in STANDBY mode.
2. Select **End Treatment** to return to the Home screen.
3. From the Home screen:
 - a. Turn the key switch to the OFF position.
 - b. Remove the key to prevent unauthorized use of the system.



NOTE


If the power cable is still connected to the electrical source, some internal circuits remain energized. To de-energize all internal circuits, turn the key switch to the OFF position and unplug the power cable from the wall.

Emergency Shut-down

If the system becomes unresponsive during laser emission, press the emergency laser stop button on the front of the console. Turn the key switch to the OFF position.


Before restarting, press the emergency laser stop button to disengage it.

Home Screen

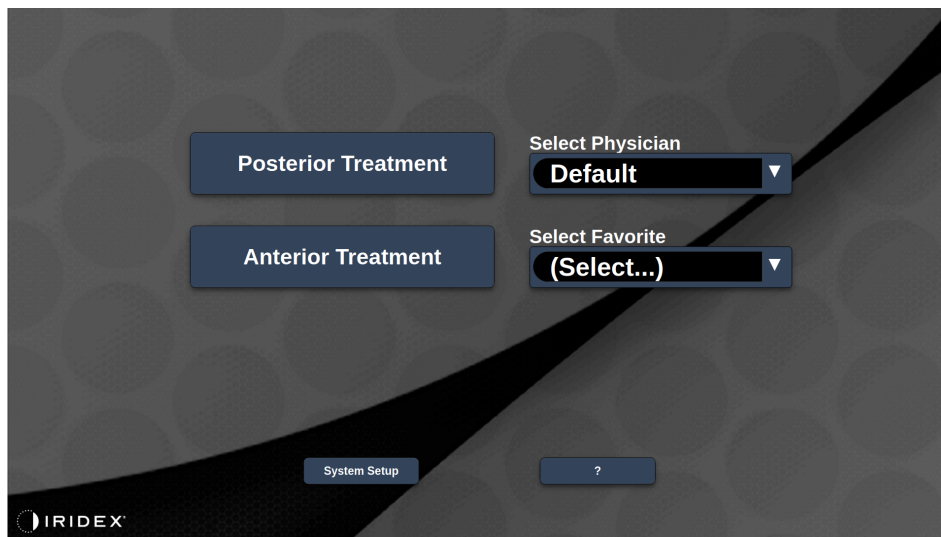


Posterior Treatment	Select to go to the Treatment screen.
Anterior Treatment	Select to go to the Treatment screen.
Select Physician	Select to display a list of physicians.
Select Favorite	Select to display a list of favorite settings.
System Setup	Select to configure the system.
?	Select to go to the system software version dialog.

After selecting the Treatment, the system will display a message at the center of the Treatment screen during system warm up.

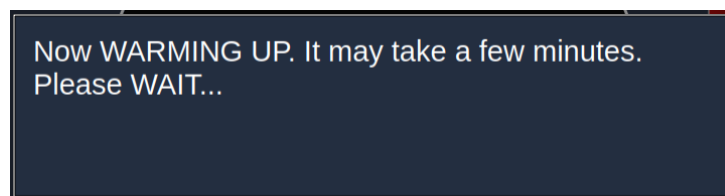

System Setup Screen

Settings Tab

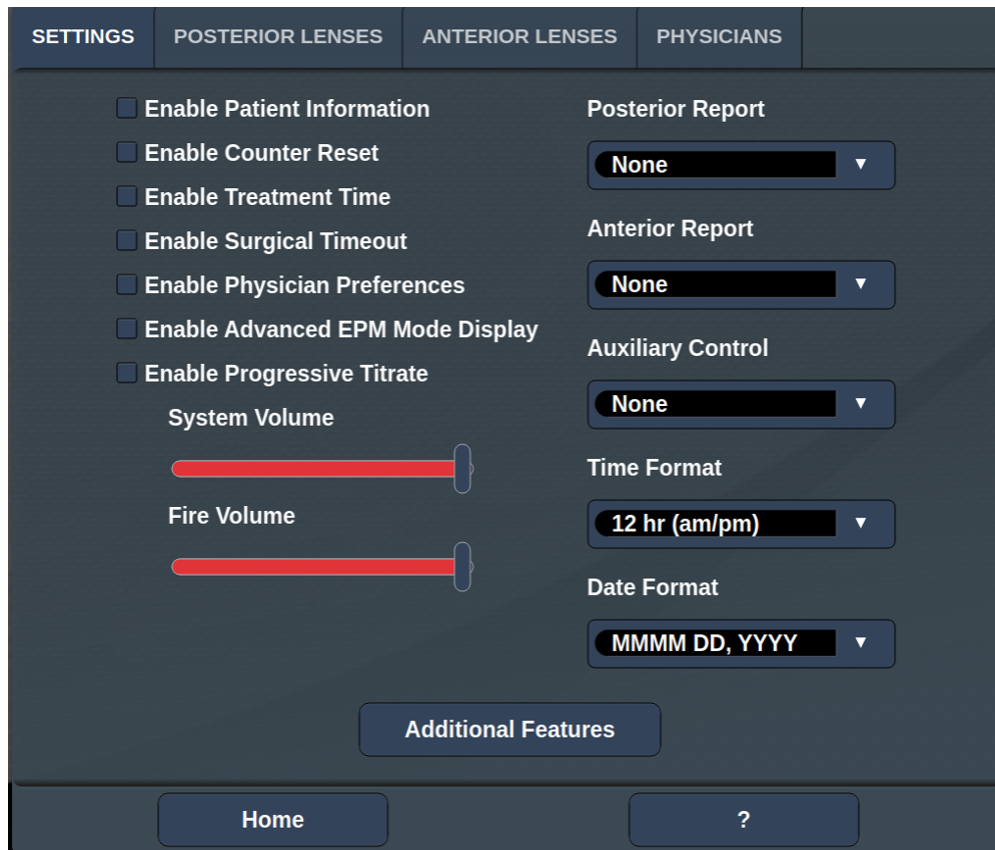


Patient Information	Enable/Disable the Patient Information pop-up prior to treatment.
Counter Reset	Enable/Disable resettable counter in the Treatment screen.
Treatment Time	Select to display the treatment time and date on the Treatment screen and report
Surgical Timeout	Enable/Disable screen to review Patient information (prior to advancing to the Treatment Screen).
Physician Preferences	Enable/Disable Physician to edit preferences.
Advanced EPM Mode Display	Enable the EPM feature in Single and Pattern Group.
Progressive Titrate	Enable/Disable Progressive titrate.
System Volume	Slide the bar left or right to adjust system volume.
Fire Volume	Slide the bar left or right to adjust fire volume.

Posterior Report	Select the export report format.
Anterior Report	Select the export report format.
Auxiliary	Select Auxiliary setting (warning light or auditory signal outside the treatment room).
Time Format	Select the time format.
Date Format	Select the date format.
Additional Features	Enter activations codes for software feature options.
Home	Select to go to the Home screen.
?	Select to go to the Help screen.

Patient Information Feature

When the Patient Information feature is enabled, the Patient Information feature allows you to include patient identification, such as patient name, number, and date of birth on the Treatment screen and in the Treatment Report.

The following window displays when you select the **Posterior Treatment** or **Anterior Treatment** button on the Home screen.

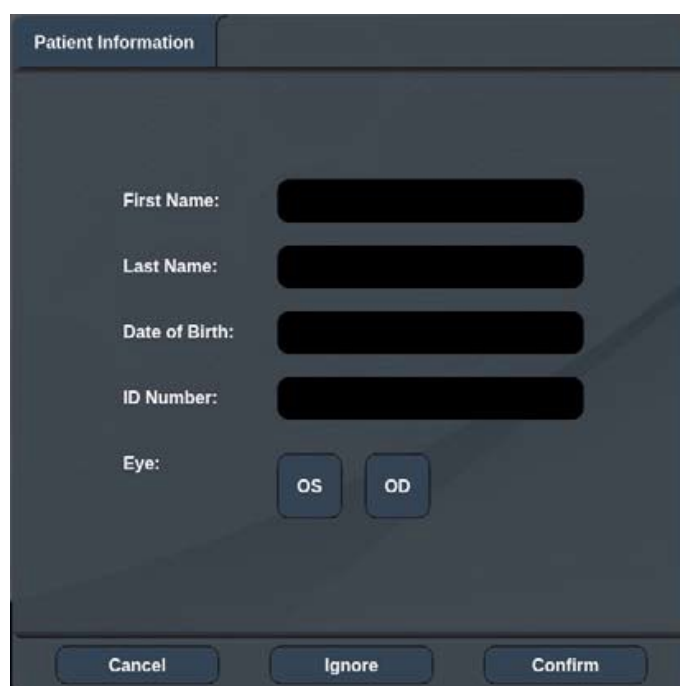

The image shows a software window titled "Patient Information". It contains several input fields: "First Name:", "Last Name:", "Date of Birth:", and "ID Number:", each followed by a dark rectangular input box. Below these is an "Eye:" label with two buttons labeled "OS" and "OD". At the bottom of the window are three buttons: "Cancel", "Ignore", and "Confirm".

The information you entered in the Patient information window is displayed on the Treatment screen and in the Treatment Report.




NOTE

*The operating system does not save this information when you select **End Treatment**.*


Counter Reset


When Counter Reset is enabled, the Treatment screen displays the Counter Control with two rows of numbers for spots and a Reset button.


- The 1st row is the treatment counter that cannot be reset.
- The 2nd row is a resettable counter.

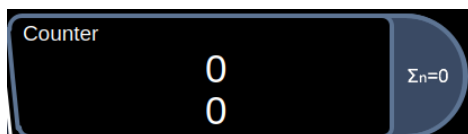

When Counter Reset is disabled, the Treatment screen displays one row of number for spots without a Reset button.

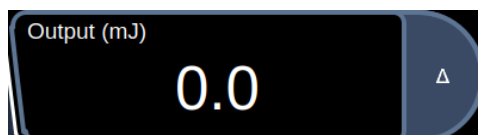


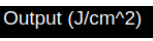If a lens is selected, selecting  will change the output in 

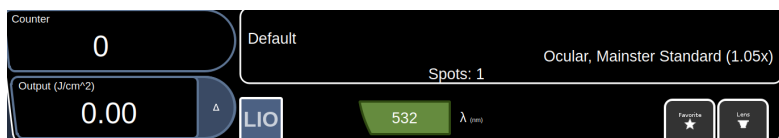

Treatment Time


When treatment time is enabled, the treatment duration is visible on the Treatment screen and in the reports.

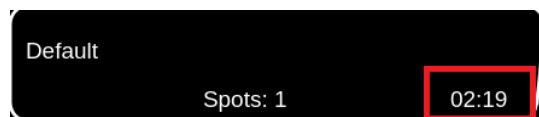

Surgical Timeout

When Enable Surgical Information is enabled, the Surgical Timeout screen is displayed prior to treatment and displays the patient information entered from the Patient Information screen (if the Patient Information feature is enabled and entered).

The surgical timeout confirmation date and time are displayed in the report when the treatment is completed.

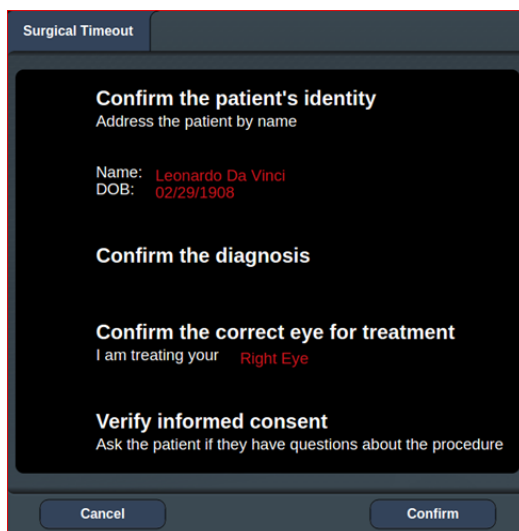


The Surgical Timeout screen is a dark-themed dialog box with a title bar labeled "Surgical Timeout". It contains four sections for confirmation:

- Confirm the patient's identity**: Address the patient by name. Name: Leonardo Da Vinci, DOB: 02/29/1908.
- Confirm the diagnosis**
- Confirm the correct eye for treatment**: I am treating your Right Eye.
- Verify informed consent**: Ask the patient if they have questions about the procedure.

At the bottom are "Cancel" and "Confirm" buttons.

Physician Preferences

When Physician Preferences is selected, Edit Preferences is enabled.

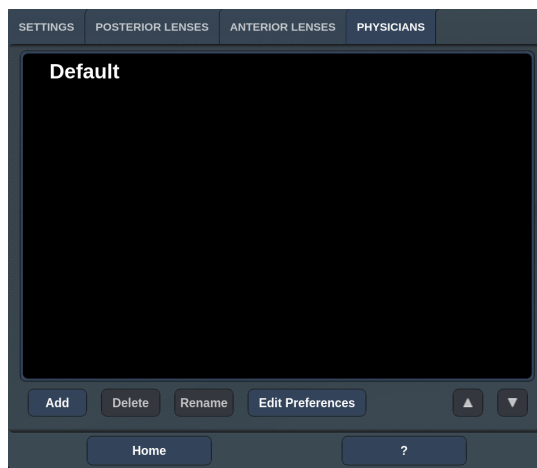

The Physician Preferences screen is a dark-themed window with a tabbed interface. The tabs are: SETTINGS, POSTERIOR LENSES, ANTERIOR LENSES, and PHYSICIANS (which is currently selected). The main area is titled "Default" and is mostly empty. At the bottom, there are buttons for "Add", "Delete", "Rename", and "Edit Preferences". Below these are "Home" and "?" buttons.

Advanced EPM Mode Display

When Advanced EPM Mode Display is enabled and has EPM Software license, EPM will be available for single and pattern group, in addition to EPM pattern group.

Progressive Titrate

When the Progressive Titrate feature is enabled, it allows the physician to output an aiming beam pattern with positions that are offset from the delivered spot with 1.5 spot diameter spacing between the titration pattern and the blinking aiming spots. If the titration pattern spots have reached the top of the FOV, the blinking aiming spots stay at the same position of the titration pattern spots. The location of titration is set to the initial position of the titration when you select the **Titrate** button to revert to the treatment mode. (See Pattern Titration under Center Panel selection for more details.)

1. Deliver the multiple spot titration pattern by activating the footswitch. Keep the footswitch activated and assess the tissue reaction with the aiming beam turned OFF.
2. When the footswitch is released, the aiming beam will shift upward by 1.5 spot size diameters and progressive shifts up each time the user decides to redo the titration.



Audio Volume

When Audio Volume is enabled, distinct audible notifications sound on an action or value change for each control activated by the user.

Fire Volume

When Fire Volume is enabled, distinct audible notifications sound when firing the laser.

Treatment Report Window

When one of the reports is selected under Posterior/Anterior report, the Treatment Report window displays when you select the **End Treatment** button at the completion of each patient treatment session.


PASCAL Posterior Laser Treatment Report

Name: _____ Patient Number: _____ DOB: _____
Date: August 17, 2022 Eye: _____ Diagnosis: _____ Procedure: _____

Patterns Used: Array
Contact Lens: None Selected

Endpoint Management [®] & Landmark [™] Settings	532 nm EM: Disabled LM: Disabled
Endpoint Management % [Median, Range]	N/A
Spot Diameter (µm)	200
Exposure (ms) [Median, Range]	20, 20 - 20
Power (mW) [Median, Range]	30, 30 - 30
Number of Spots [Total, (#LM, #EM)]	17 (17, 0)
Energy (mJ) [Median, Range]	5, 4 - 5
Fluence (J/cm²) [Median, Range]	No Lens Selected

Notes:

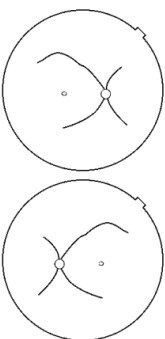


1 of 1
100 %
EXPORT
CLOSE

	Select right arrow or left arrow  to show the next or previous page.
	Select  to zoom in/out in the preview.
Export	Select to export the Treatment Report to the connected USB drive.
Close	Select to close the Treatment Report window and go to the Home screen.



NOTE

- When you select the **CLOSE** button, the system returns to the Home screen and the current session treatment report is no longer available.
- The Export is disabled if no USB drive is plugged into the table's USB port. If the connected USB drive doesn't meet the requirements, the **EXPORT** button is still disabled, and the following message pops up.

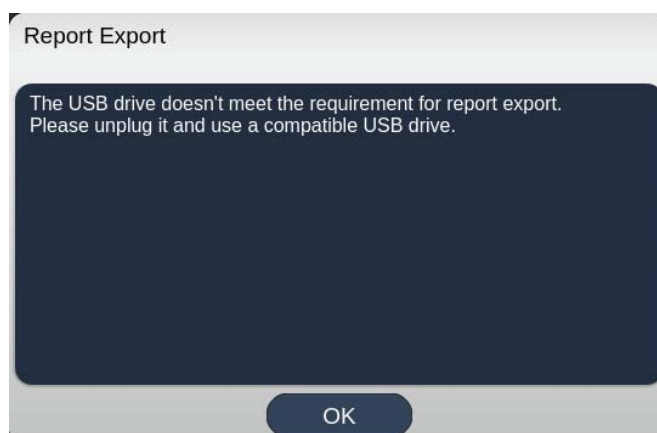


All the reports exported from Iridex Pascal will be stored in the 'synthesis_report' folder in the USB drive. You can plug the USB drive into a PC to copy or delete the exported reports.

The USB drive should be dedicated for the Report Export purpose. Any other usage will make it incompatible for Report Export, in which case you will have to reformat it according to the [Procedure to prepare the USB drive before the first time use](#).

Procedure to prepare the USB drive before the first time use

1. Purchase a new USB drive of USB 2.0 or USB 3.0 with at least 8GB in size.
2. On a Windows 10 PC, open Windows Explorer, plug in the USB drive, and identify the newly added drive. Rename the drive to "SYNTHESIS"

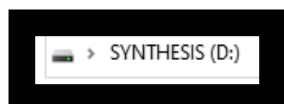

Posterior Reports

Select the report setting from the following: None, Basic, or Advanced

PASCAL Posterior Laser Treatment Report


Name: _____ Patient Number: _____ DOB: _____
 Date: _____ Eye: _____ Diagnosis: _____ Procedure: _____

Patterns Used: Triple Arc, A+B Octants (15 ms), Hexagon (15 ms), Octants
 Contact Lens: Volk, SuperQuad 160 (2,00x), Ocular, Fundus Laser (1,08x), Rodenstock Schlegel
 Panfunduscope (1,50x)

577 nm
Endpoint Management (EM) & Landmark (LM) Settings
EM: Disabled LM: Disabled
Endpoint Management % [Median, Range]
N/A
Spot Diameter (µm)
100, 200
Exposure (ms) [Median, Range]
15, 10 - 30
Power (mW) [Median, Range]
100, 0 - 1000
Number of Spots [Total, (#LM, #EM)]
771 (771, 0)
Energy (mJ) [Median, Range]
18, 0 - 960
Fluence (J/cm²) [Median, Range]
4.1, 0.0 - 81.9

Treatment Time: 30:00 (mm:ss)
 Note:

_____/_____/_____ MD
 _____ Signature
 _____ / _____ / _____ Date

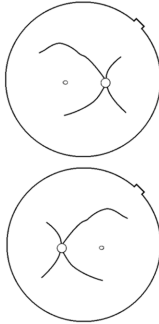

Sample Posterior Basic Treatment Report

PASCAL Posterior Laser Treatment Report


Name: _____ Patient Number: _____ DOB: _____
 Date: _____ Eye: _____ Diagnosis: _____ Procedure: _____

Total # of Spots: 784 Patterns Used: Triple Arc, arc, A+B Octants, Octants (15 ms), Hexagon (15 ms), Octants
 Power (mW) 100 Spot Diameter 100, 200 Exposure (ms) 15,
 [Median, Range]: 0 - 1000 (µm); [Median, Range]: 10 - 30
 Contact Lens: Volk, SuperQuad 160 (2,00x), Ocular, Fundus Laser (1,08x), Rodenstock Schlegel
 Panfunduscope (1,50x)

A	Patterns	# Spots	Power (mW) [Median, Range]	Exposure (ms) [Median, Range]	EM % [Median, Range]	Spot Spacing (°) [Median, Range]	Energy (mJ) [Median, Range]
S77 nm	Triple Arc	128 (128, 0)	300, 225 - 300	25, 20 - 20	N/A	0.5, 1.5	10, 4 - 196
	arc	100 (100, 0)	1000, 1000 - 1000	30, 30 - 30	N/A	0.25	950, 810 - 1200
	A+B Octants	112 (112, 0)	100, 100 - 100	15, 15 - 15	N/A	1.5	14, 14 - 14
	Octants (15ms)	208 (16, 160)	100, 100 - 100	15, 15 - 15	30, 30 - 30	0.25	14, 14 - 14
	Hexagon (15 ms)	85 (80, 85)	100, 100 - 100	15, 15 - 15	30, 30 - 30	2	15, 15 - 15
	Octants	144 (144, 0)	0, 0 - 0	15, 15 - 15	N/A	1.5	0, 0 - 0

Total Area Treated (mm²): 31.38 Total Energy Delivered (mJ): 4004.5
 Treatment Time: 00:30 minutes Fluence (J/cm²) 1.2, 0.0 - 81.9
 [Median, Range]:

Endpoint Management (EM) & Landmark (LM) Settings EM = On, LM = On
 _____ MD
 _____ Signature
 _____ / _____ / _____ Date

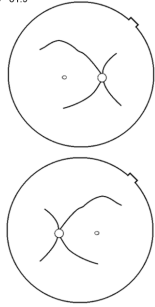

Sample Posterior Advanced Treatment Report

Anterior Reports

Select the report setting from the following: None or Basic

PASCAL Anterior Laser Treatment Report

Name: _____ Patient Number: _____ DOB: _____
 Date: _____ Eye: _____ Diagnosis: _____ Procedure: _____

Patterns Used: PSLT, Array (Ant.)
 Contact Lens: Ocular, Latina SLT Gonio (1.00x)


	577 nm Titration	577 nm Summary
Spot Diameter (µm)	100	100
Exposure (ms) [Median, Range]	10, 10 - 10	5, 5 - 20
Power (mW) [Median, Range]	85, 70 - 100	100, 100 - 200
Number of Spots	4	1278
Energy (mJ) [Median, Range]	1, 1 - 1	20, 6 - 20
Fluence (J/cm²) [Median, Range]	10.8, 8.9 - 12.7	6.4, 6.4 - 25.5
Treatment Plan		360°
Total Angle Treated		360°

Note:

_____ MD

_____ Signature

____ / ____ / ____ Date

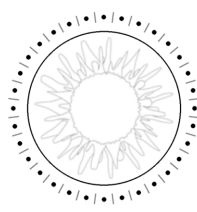

Sample Anterior Basic Treatment Report



NOTE

The graphical depiction of the area treated on the Anterior Laser Treatment Report represents the physical area of the Trabecular Meshwork (TM) that was treated during the procedure. If any areas of the TM were retreated during a treatment, these areas will be visually depicted with a solid circle.

Auxiliary Control

Aux setting to enable Aux

Auxiliary Control None On Ready On Treat	None	Aux off all the time
	Ready	Lid up Aux when changed from STANDBY to READY mode.
	Treat	Lid up Aux when changed from STANDBY to TREAT.

Time Format

Select the time format for the time of day displayed in current local time.

Date Format

Select the date format for the date of birth and treatment date on the report, and date of birth on the Surgical Timeout confirmation screen.


Additional Features

Enter Iridex Pascal License Code to enable the license key of MicroPuse/EPM/PSLT.

Lens Tab

Posterior/Anterior Lenses

Each type of treatment consists of a List box that displays the Contact List. Each Contact Lens List & Custom Contact Lens Magnification are independent from the Posterior's/Anterior's.

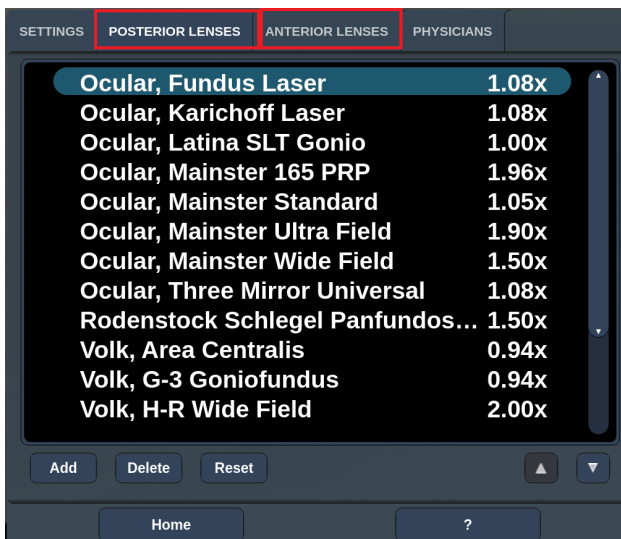

Custom Contact Lens Magnification

Select **Add** on the Lens tab and type the information in the Enter lens Information pop-up dialog.

Physicians Tab

The Physicians tab allows you to add, delete, or reset a set list of physicians. When Physician Preference is enabled in settings, Physician Preference can be renamed or edited.

Physician Preferences: tester

Ocular Settings

☐ Ocular Setting Reminder

OS OD

0 0

Advance Direction

Clockwise

Auto-Transition after Enhanced Octants

Array

Cancel Ok

Ocular Setting Reminder	Displays settings to remind the user to set ocular power before treatment.
Advance Direction	Select direction for the Auto Advance feature.
Auto-Transition after Enhanced Octants	Will advance to pattern after Enhanced Octant is completed (available when Endpoint Management is activated).

Auto-Transition Feature

If Array or Hexagon pattern is selected in Auto-Transition after Enhanced Octants in the Physician Preferences screen, it allows the user to transition to the selected pattern automatically after treating the last subset of EpM Enhanced Octants and to STANDBY mode.



In the case of Array pattern selected




In the case of Hexagon pattern selected

Auto-Transition After Enhanced Octants

Treatment Screen

Top Panel

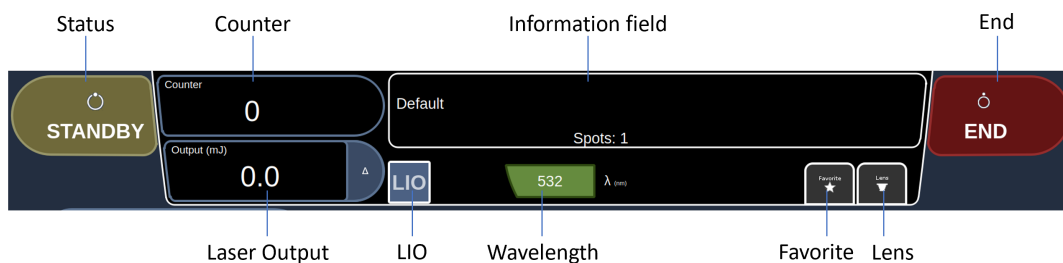

Status	Select the system status (STANDBY or READY mode).
Counter	Displays the number of applications during treatment.
	Σn=0 - Select to reset the counter (Enable in preferences).
Information	Displays optional information and warning messages.
End	Select to end treatment.
Output	Displays laser output value.
	Δ - Select to toggle units.
LIO	Select LIO.
Wavelength	532 or 577 as fixed value, depending on the system.
Favorite	Select to open the Favorite window. See " Favorite Window ".
Lens	Select to open the Contact Lens menu.


Favorite Window

Select the favorite icon in the Posterior/Anterior Treatment screen to open the Favorite window.


Display favorite details information (pattern, power, exposure, spot diameter, spacing, lens) when the favorite entry is pressed and held for 1 second in the Favorite list. To close the pop-up, touch the pop-up or anywhere on the screen. If the favorite is pressed and released shorter than 1 second, the favorite will be loaded.



Select the **Add** button and enter the name for the favorite to save the current parameter setting as your favorite.

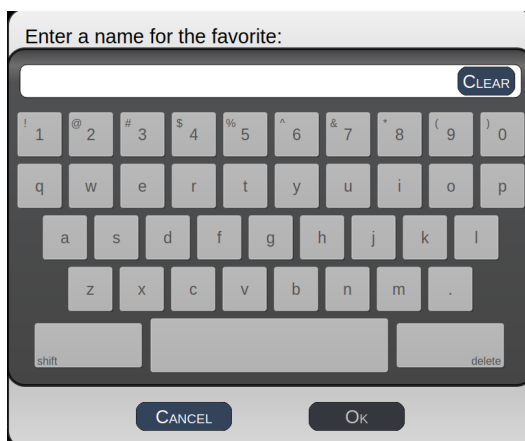

When a favorite is loaded, the favorite name will be displayed in the status area.

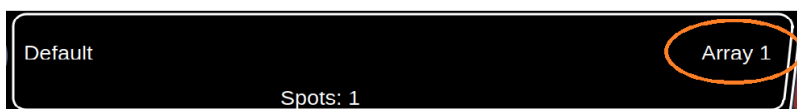

NOTE

If any parameter is changed, the favorite name in the status area will be cleared.

Edit Favorites Window

Select the **Edit** button to enter Edit mode.




When a favorite is selected, the **Rename**, **Remove**, **Save**, and **Move** buttons are enabled.

Done	Select when you are finished editing.
Rename	Select to change the name of the selected favorite.
Remove	Select to remove the selected favorite from Favorites window.
Save	Select to save the edit for the selected favorite.


Move	Select to reorder the favorite list.
Close	Select to close the Favorite window.

To move the location of a favorite list:

1. Select **Edit**.
2. Select the favorite you want to move.
3. Select **Move**.
4. Drag the list to the location you want to move to.

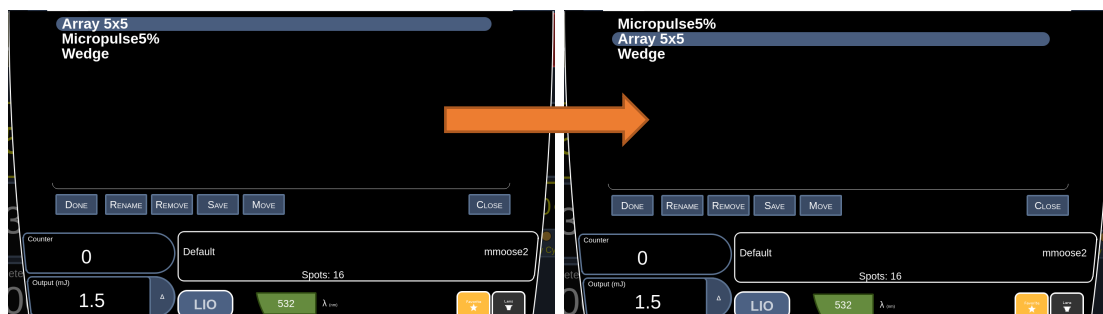


Center Panel

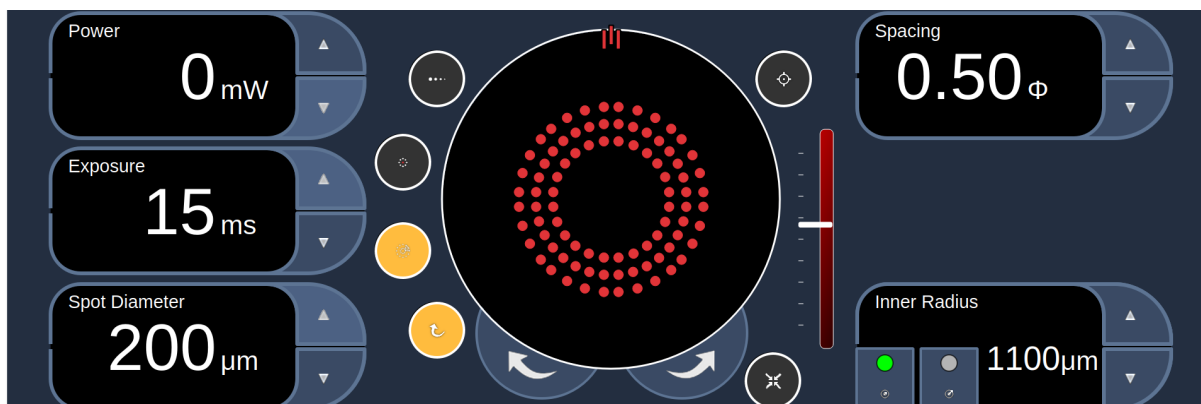

Controls

Power

Description

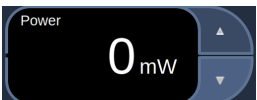
	Up/Down arrow button to adjust Power, or turn the Power Knob on the slit lamp.
---	--

Range

Spot Diameter	Range
LIO	Power Range: 30 mW - 1500 mW
50 µm	Power Range: 30 mW - 1000 mW
100 µm	Power Range: 30 mW - 1500 mW
200 µm	Power Range: 30 mW - 2000 mW
300 µm	Power Range: 30 mW - 1500 mW

Exposure

Description

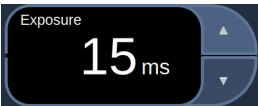	Up/Down arrow button to adjust Duration of the treatment laser delivered to the target.
---	---

Range: Posterior

Pattern (Default)		Range
Single Spot	LIO	Exposure Range: 10 ms - 1000 ms
	Single Spot (10 ms)	Exposure Range: 10 ms - 1000 ms
Pattern	2x2 (20 ms)	Exposure Range: 10 ms - 80 ms
	3x3 (20 ms)	Exposure Range: 15 ms - 80 ms
	4x4, 5x5 (20 ms)	Exposure Range: 15 ms - 30 ms
	If n <= 2 and m <= 2 Min: 10 ms, If n > 2 or m > 2 Min: 15 ms, If n <= 3 and m <= 3 Max: 80 ms, If n > 3 or m > 3 Max: 30 ms	
	Hexagon (20 ms)	Exposure Range: 15 ms - 30 ms
	Triple Arc (20 ms)	Exposure Range: 15ms - 30 ms
	Single Spot (20 ms)	Exposure Range: 10 ms - 1000 ms
	Wedge (20 ms)	Exposure Range: 15 ms - 30 ms
	Arc (20 ms)	Exposure Range: 15 ms - 30 ms
	Triple Ring (20 ms)	Exposure Range: 15 ms – 20 ms
	Line (20 ms)	Exposure Range: 15 ms - 30 ms
	Octant (10 ms)	Exposure Range: 10 ms - 15 ms
	Enhanced Octant (10 ms)	Exposure Range: 10 ms - 15 ms


Pattern (Default)		Range	
MicroPulse	2x2 (20 ms)	Exposure Range:	10 ms - 300 ms
	3x3 (20 ms)	Exposure Range:	15 ms - 300 ms
	4x4, 5x5 (20 ms)	Exposure Range:	15 ms - 300 ms
	Hexagon (20 ms)	Exposure Range:	15 ms - 300 ms
	Triple Arc (20 ms)	Exposure Range:	15 ms - 300 ms
	Single Spot (20 ms)	Exposure Range:	10 ms - 1000 ms
	Wedge (20 ms)	Exposure Range:	15 ms - 300 ms
	Arc (20 ms)	Exposure Range:	15 ms - 300 ms
	Triple Ring (20 ms)	Exposure Range:	15 ms - 300 ms
	Line (20 ms)	Exposure Range:	15 ms - 300 ms
	Octant (10 ms)	Exposure Range:	10 ms - 300 ms
	Enhanced Octant (10 ms)	Exposure Range:	10 ms - 300 ms
EpM	Enhanced Octant (15 ms)	Exposure Range:	10 ms - 15 ms
	Array 2x2 (15 ms)	Exposure Range:	10 ms - 80 ms
	Array 3x3 (15 ms)	Exposure Range:	15 ms - 80 ms
	Array 4x4 (15 ms)	Exposure Range:	15 ms - 30 ms
	If n <= 2 and m <= 2 Min: 10 ms, If n > 2 or m > 2 Min: 15 ms, If n <= 3 and m <= 3 Max: 80 ms, If n > 3 or m > 3 Max: 30 ms		
	Hexagon (15 ms)	Exposure Range:	15 ms - 30 ms
	Single Spot (15 ms)	Exposure Range:	10 ms - 80 ms

Range Anterior

Pattern		Range	
Anterior	PSLT (5 ms)	Exposure Range:	5ms
	Array 2x2 (20 ms)	Exposure Range:	10 ms - 80 ms
	Array 3x3 (20 ms)	Exposure Range:	15 ms - 80 ms
	If n <= 2 and m <= 2 Min: 10 ms, If n > 2 or m > 2 Min: 15 ms, If n <= 3 and m <= 3 Max: 80 ms, If n > 3 or m > 3 Max: 30 ms		
	LIO	Exposure Range:	10 ms - 1000 ms
	Single Spot (10 ms)	Exposure Max:	10 ms - 1000 ms

Spot Diameter Control

Description

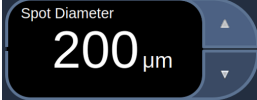	Up/Down arrow button to adjust Spot size of the treatment laser delivered to the target
---	---

Range

Pattern (yellow/green laser)		Spot Diameter available (Y/N)?			
		50 μm	100 μm	200 μm	300 μm
Single Spot	LIO	N	N	N	N
	Single Spot	Y	Y (Default)	Y	Y
Pattern	Array 2X2	Y	Y	Y (Default)	Y
	Array 3X3, 4X4, 5X5	N	Y	Y (Default)	Y
	Hexagon	N	Y	Y (Default)	Y
	Triple Arc	N	N	Y (Default)	Y
	Single Spot	Y	Y (Default)	Y	Y
	Wedge	N	Y	Y (Default)	Y
	Arc	N	N	Y (Default)	Y
	Triple Ring	N	N	Y (Default)	Y
	Line	N	Y	Y (Default)	Y
	Octant	N	Y (Default)	Y	N
	Enhanced Octant	N	Y (Default)	Y	N
MicroPulse	The same as Pattern Pallet				
EpM	Enhanced Octant	N	Y	Y (Default)	N
	Array 2x2	Y	Y	Y (Default)	Y
	Array 3x3, 4x4	N	Y	Y (Default)	Y
	Hexagon	N	Y	Y (Default)	Y
	Single Spot	N	Y (Default)	Y	N

Spacing Field

Description

	Up/Down arrow button to adjust the edge to edge spacing of spots in the output treatment or aim laser pattern
---	---

Range


Pattern Pallet Array (Default: 0.50Φ)	Spot Diameter	Range
2x2	50 μm	Spacing Range: 0.00Φ - 3.00Φ
2x2, 3x3, 4x4, 5x5	100 μm	Spacing Range: 0.00Φ - 3.00Φ
2x2, 3x3, 4x4, 5x5	200 μm	Spacing Range: 0.00Φ - 3.00Φ
2x2, 3x3, 4x4, 5x5	300 μm	Spacing Range: 0.00Φ - 1.50Φ
Pattern Pallet Hexagon (Default: 0.50Φ)	Spot Diameter	Range
7 Spot, 19 Spot, 37 Spot	100 μm	Spacing Range: 0.00Φ - 3.00Φ
7 Spot, 19 Spot	200 μm	Spacing Range: 0.00Φ - 3.00Φ
37 Spot		Spacing Range: 0.00Φ - 2.50Φ
7 Spot	300 μm	Spacing Range: 0.00Φ - 3.00Φ
19 Spot		Spacing Range: 0.00Φ - 1.50Φ
37 Spot		Spacing Range: 0.00Φ - 0.50Φ
Pattern Pallet Triple Arc (Default: 0.25Φ)	Spot Diameter	Range
(0.25Φ)	200 μm	Spacing Range: 0.25Φ - 3.00Φ
	300 μm	Spacing Range: 0.25Φ - 1.50Φ
Pattern Pallet Wedge (Default: 0.50Φ)	Spot Diameter	Range
6,10,15,21 Spot	100 μm	Spacing Range: 0.00Φ - 3.00Φ
6,10 Spot	200 μm	Spacing Range: 0.00Φ - 3.00Φ
15 Spot		Spacing Range: 0.00Φ - 2.50Φ
21 Spot		Spacing Range: 0.00Φ - 2.00Φ
6,10 Spot	300 μm	Spacing Range: 0.00Φ - 1.00Φ
15 Spot		Spacing Range: 0.00Φ - 0.75Φ
21 Spot		Spacing Range: 0.00Φ - 0.50Φ
Pattern Pallet Triple Ring (Default: 0.50Φ)	Spot Diameter	Range
Radius: 500μm	200 μm	Spacing Range: 0.50Φ
Radius: 500μm - 600μm	300 μm	Spacing Range: 0.00Φ - 0.50Φ
Radius: 700μm		Spacing Range: 0.25Φ - 0.50Φ
Radius: 800μm - 1000μm		Spacing Range: 0.50Φ
Pattern Pallet Arc (Default: 0.50Φ)	Spot Diameter	Range

Pattern Pallet Array (Default: 0.50Φ)		Spot Diameter	Range
Radius: 1300 μm – 1500 μm Radius: 1100 μm – 1200 μm Radius: 500 μm – 1000 μm		200 μm	Spacing Range: 0.50Φ Spacing Range: 0.25Φ - 0.50Φ Spacing Range: 0.00Φ - 0.50Φ
Radius: 500 μm – 2000 μm		300 μm	Spacing Range: 0.00Φ - 0.50Φ
Pattern Pallet Line (Default: 0.50Φ)		Spot Diameter	Range
2 - 10 spots		100 μm	Spacing Range: 0.00Φ - 3.00Φ
10 spots 2 - 9 spots		200 μm	Spacing Range: 0.00Φ - 1.25Φ Spacing Range: 0.00Φ - 1.50Φ
10 spots 9-8 spots 7 spots 6 spots 5 spots 4 spots 3 spots		300 μm	Spacing Range: 0.00Φ Spacing Range: 0.00Φ - 0.25Φ Spacing Range: 0.00Φ - 0.50Φ Spacing Range: 0.00Φ - 1.00Φ Spacing Range: 0.00Φ - 1.50Φ Spacing Range: 0.00Φ - 2.00Φ Spacing Range: 0.00Φ - 3.00Φ
Pattern Pallet Octant (Default: 1.50Φ)		Spot Diameter	Range
Radius: 1100 μm (# Octants: 1) Radius: 1200 μm (# Octants: 1) Radius: 1300 μm - 1400μm (# Octants: 1) Radius: 1500 μm (# Octants: 1) Radius: 1600 μm (# Octants: 1)		100 μm	Spacing Range: 0.25Φ - 2.00Φ Spacing Range: 0.25Φ - 1.50Φ Spacing Range: 0.25Φ - 1.00Φ Spacing Range: 0.25Φ - 0.50Φ Spacing Range: 0.25Φ
Radius: 1100 μm (# Octants: 1) Radius: 1200 μm (# Octants: 1)		200 μm	Spacing Range: 0.25Φ - 0.50Φ Spacing Range: 0.25Φ
Pattern Pallet Enhanced Octant (Default: 1.50Φ)		Spot Diameter	Range
Outer Radius:	Inner Radius Range:		
2200 μm	500 μm - 1700 μm	100 μm	Spacing Range: 0.00Φ - 3.00Φ
2200 μm	500 μm - 1800 μm		Spacing Range: 0.00Φ - 2.00Φ
2200 μm	500 μm - 1900 μm		Spacing Range: 0.00Φ - 1.00Φ
2200 μm	500 μm - 2000 μm		Spacing Range: 0.00Φ - 0.00Φ
2200 μm	500 μm - 1600 μm	200 μm	Spacing Range: 0.00Φ - 1.00Φ
2200 μm	500 μm - 1700 μm		Spacing Range: 0.00Φ - 0.50Φ
2200 μm	500 μm - 1800 μm		Spacing Range: 0.00Φ - 0.00Φ
MicroPulse Pallet Array		Spot Diameter	Range
The same as Pattern Pallet			
EpM Pallet Enhanced Octant (Default: 0.25Φ)		Spot Diameter	Range
Outer Radius:	Inner Radius Range:		
2200 μm	500 μm - 1700 μm	100 μm	Spacing Range: 0.00Φ - 3.00Φ
2200 μm	500 μm - 1800 μm		Spacing Range: 0.00Φ - 2.00Φ

Pattern Pallet Array (Default: 0.50Φ)		Spot Diameter	Range
2200 μm	500 μm - 1900 μm		Spacing Range: 0.00Φ - 1.00Φ
2200 μm	500 μm - 2000 μm		Spacing Range: 0.00Φ - 0.00Φ
2200 μm	500 μm - 1600 μm	200 μm	Spacing Range: 0.00Φ - 1.00Φ
2200 μm	500 μm - 1700 μm		Spacing Range: 0.00Φ - 0.50Φ
2200 μm	500 μm - 1800 μm		Spacing Range: 0.00Φ - 0.00Φ
EpM Pallet Array (Default: 0.25Φ)		Spot Diameter	Range
2x2, 3x3, 4x4		50 μm	Spacing Range: 0.00Φ - 3.00Φ
2x2, 3x3, 4x4		100 μm	Spacing Range: 0.00Φ - 3.00Φ
2x2, 3x3, 4x4		200 μm	Spacing Range: 0.00Φ - 3.00Φ
2x2, 3x3, 4x4		300 μm	Spacing Range: 0.00Φ - 1.50Φ
EpM Pallet Hexagon (Default: 0.25Φ)		Spot Diameter	Range
7 Spot, 19 Spot, 37 Spot		100 μm	Spacing Range: 0.00Φ - 3.00Φ
7 Spot, 19 Spot		200 μm	Spacing Range: 0.00Φ - 3.00Φ
37 Spot			Spacing Range: 0.00Φ - 2.50Φ
7 Spot		300 μm	Spacing Range: 0.00Φ - 3.00Φ
19 Spot			Spacing Range: 0.00Φ - 1.50Φ
37 Spot			Spacing Range: 0.00Φ - 0.50Φ
ANTERIOR Pallet Array (Default: 0.50Φ)		Spot Diameter	Range
2x2		50 μm	Spacing Range: 0.00Φ - 3.00Φ
2x2, 3x3		100 μm	Spacing Range: 0.00Φ - 3.00Φ
2x2, 3x3		200 μm	Spacing Range: 0.00Φ - 3.00Φ
2x2, 3x3		300 μm	Spacing Range: 0.00Φ - 1.50Φ

Endpoint Management

Description

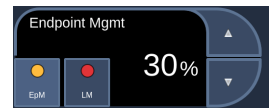	Select EpM to change the value of percentage range
---	--

Range

Pattern (Default)	Range
EpM	10% to 95%.

MicroPulse

Description


	• Select Duration to manually set the value. • Select Interval to manually set the value. • Select Duty Cycle to toggle between default values of 5%,10%,15%.
---	---

Range

Pattern (Default)	Range
Duration	Exposure Range: 0.05 – 1.00 ms
Interval	Exposure Range: 0.10 – 10.00 ms

Radius Field

Description

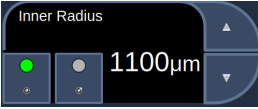	• Up/Down arrow button to increment or decrement radius • Left side Radio button to adjust Inner Radius • Right side Radio button to adjust Outer Radius
--	--


Range

Spot Spacing	Spot Diameter	Range
Pattern Pallet Triple Ring (Default Radius: 500 µm)		
Spacing: 0.50Φ	200 µm	Radius Range: 500 µm
Spacing: 0.00Φ	300 µm	Radius Range: 500 µm - 600 µm
Spacing: 0.25Φ		Radius Range: 500 µm - 700 µm
Spacing: 0.50Φ		Radius Range: 500µm - 1000 µm
Pattern Pallet Arc (Default Radius: 1000 µm)		
Spacing: 0.00Φ	200 µm (Auto Advance On)	Radius Range: 500 µm - 1000 µm
Spacing: 0.25Φ		Radius Range: 500 µm - 1200 µm
Spacing: 0.50Φ		Radius Range: 500 µm - 1500 µm
Spacing: 0.00Φ	200 µm (Auto Advance Off)	Radius Range: 500 µm - 1400 µm
Spacing: 0.25Φ		Radius Range: 500 µm - 1600 µm
Spacing: 0.50Φ		Radius Range: 500 µm - 1900 µm
Spacing: 0.00Φ - 0.50 Φ	300 µm	Radius Range: 500 µm - 2000 µm
Pattern Pallet Octant (Default Radius: 1100 µm)		
Spacing: 0.25Φ	100 µm	Radius Range: 1100 µm - 1600 µm

Spot Spacing	Spot Diameter	Range
Spacing: 0.50Φ		Radius Range: 1100 μm - 1500 μm
Spacing: 1.00Φ		Radius Range: 1100 μm - 1400 μm
Spacing: 1.50Φ		Radius Range: 1100 μm - 1200 μm
Spacing: 2.00Φ		Radius Range: 1100 μm
Spacing: 0.25Φ	200 μm	Radius Range: 1100 μm - 1200 μm
Spacing: 0.50Φ		Radius Range: 1100 μm
Pattern Pallet Enhanced Octant (Default Outer Radius: 2000 μm)		
Spacing: 0.00Φ	100 μm Set Outer Radius to 2200 μm	Inner Radius Range: 500 μm - 2000 μm Outer Radius Max: 2200 μm
Spacing: 0.25Φ-1.00Φ		Inner Radius Range: 500 μm - 1900 μm Outer Radius Max: 2200 μm
Spacing: 1.50Φ-2.00Φ		Inner Radius Range: 500 μm - 1800 μm Outer Radius Max: 2200 μm
Spacing: 2.50Φ-3.00Φ		Inner Radius Range: 500 μm - 1700 μm Outer Radius Max: 2200 μm
Spacing: 0.00Φ	200 μm Set Outer Radius to 2200 μm	Inner Radius Range: 500 μm - 1800 μm Outer Radius Max: 2200 μm
Spacing: 0.25Φ-0.50Φ		Inner Radius Range: 500 μm - 1700 μm Outer Radius Max: 2200 μm
Spacing: 1.00Φ		Inner Radius Range: 500 μm - 1600 μm Outer Radius Max: 2200 μm
MicroPulse Pallet Triple Ring (Default Radius: 500 μm)		
The same as Pattern Pallet		
EpM Enhanced Octant (Default Outer Radius: 2200 μm)		
Spacing: 0.00Φ	100 μm Set Outer Radius to 2200 μm	Inner Radius Range: 500 μm - 2000 μm Outer Radius Max: 2200 μm
Spacing: 0.25Φ-1.00Φ		Inner Radius Range: 500 μm - 1900 μm Outer Radius Max: 2200 μm
Spacing: 1.50Φ-2.00Φ		Inner Radius Range: 500 μm - 1800 μm Outer Radius Max: 2200 μm
Spacing: 2.50Φ-3.00Φ		Inner Radius Range: 500 μm - 1700 μm Outer Radius Max: 2200 μm
Spacing: 0.00Φ	200 μm Set Outer Radius to 2200 μm	Inner Radius Range: 500 μm - 1800 μm Outer Radius Max: 2200 μm
Spacing: 0.25Φ-0.50Φ		Inner Radius Range: 500 μm - 1700 μm Outer Radius Max: 2200 μm
Spacing: 1.00Φ		Inner Radius Range: 500 μm - 1600 μm Outer Radius Max: 2200 μm

Curvature Field

Description

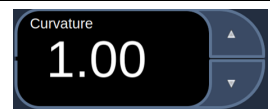
	Up/Down arrow button to adjust the curvature
---	--

Range

Pattern	Range
Triple Arc	0.00, 0.25, 0.50, 0.75, 1.00, 1.25, 1.50, 1.75, 2.00
PSLT	0.00, 0.25, 0.50, 0.75, 1.00, 1.25, 1.50, 1.75, 2.00, 2.50, 3.00

Rep Rate Field

Description

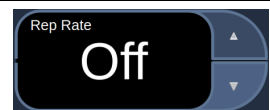	Up/Down arrow button to adjust repetition rate
---	--

Range

Pattern	Range
Single Spot in Single Spot Pallet	Off – Max Exposure is 1000 ms (Default) 1.0Hz – Max Exposure is 800 ms 1.5Hz – Max Exposure is 500 ms 2.0Hz – Max Exposure is 400 ms 3.0Hz – Max Exposure is 250 ms 4.0Hz – Max Exposure is 200 ms 5.0Hz – Max Exposure is 150 ms 6.0Hz – Max Exposure is 100 ms 7.0Hz – Max Exposure is 100 ms 8.0Hz – Max Exposure is 100 ms
Single Spot in Pattern Pallet	Off – Max Exposure is 1000 ms (Default) 1.0Hz – Max Exposure is 800 ms 1.5Hz – Max Exposure is 500 ms 2.0Hz – Max Exposure is 400 ms 3.0Hz – Max Exposure is 250 ms 4.0Hz – Max Exposure is 200 ms 5.0Hz – Max Exposure is 150 ms 6.0Hz – Max Exposure is 100 ms 7.0Hz – Max Exposure is 100 ms 8.0Hz – Max Exposure is 100 ms

Pattern	Range
Single Spot in MicroPulse Pallet	Off – Max Exposure is 1000 ms (Default) 1.0Hz – Max Exposure is 800 ms 1.5Hz – Max Exposure is 500 ms 2.0Hz – Max Exposure is 400 ms 3.0Hz – Max Exposure is 250 ms 4.0Hz – Max Exposure is 200 ms 5.0Hz – Max Exposure is 150 ms 6.0Hz – Max Exposure is 100 ms 7.0Hz – Max Exposure is 100 ms 8.0Hz – Max Exposure is 100 ms
Single Spot in EpM Pallet	No Repetition Rate
Single Spot in ANT	Off – Max Exposure is 1000 ms (Default) 1.0Hz – Max Exposure is 800 ms 1.5Hz – Max Exposure is 500 ms 2.0Hz – Max Exposure is 400 ms 3.0Hz – Max Exposure is 250 ms 4.0Hz – Max Exposure is 200 ms 5.0Hz – Max Exposure is 150 ms 6.0Hz – Max Exposure is 100 ms 7.0Hz – Max Exposure is 100 ms 8.0Hz – Max Exposure is 100 ms
LIO in Single Spot Pallet	Off – Max Exposure is 1000 ms (Default) 1.0Hz – Max Exposure is 800 ms 1.5Hz – Max Exposure is 500 ms 2.0Hz – Max Exposure is 400 ms 3.0Hz – Max Exposure is 250 ms 4.0Hz – Max Exposure is 200 ms 5.0Hz – Max Exposure is 150 ms 6.0Hz – Max Exposure is 100 ms 7.0Hz – Max Exposure is 100 ms 8.0Hz – Max Exposure is 100 ms
LIO in MicroPulse Pallet	Off – Max Exposure is 1000 ms (Default) 1.0Hz – Max Exposure is 800 ms 1.5Hz – Max Exposure is 500 ms 2.0Hz – Max Exposure is 400 ms 3.0Hz – Max Exposure is 250 ms 4.0Hz – Max Exposure is 200 ms 5.0Hz – Max Exposure is 150 ms 6.0Hz – Max Exposure is 100 ms 7.0Hz – Max Exposure is 100 ms 8.0Hz – Max Exposure is 100 ms
LIO in ANT	Off – Max Exposure is 1000 ms (Default) 1.0Hz – Max Exposure is 800 ms 1.5Hz – Max Exposure is 500 ms 2.0Hz – Max Exposure is 400 ms 3.0Hz – Max Exposure is 250 ms 4.0Hz – Max Exposure is 200 ms 5.0Hz – Max Exposure is 150 ms 6.0Hz – Max Exposure is 100 ms 7.0Hz – Max Exposure is 100 ms 8.0Hz – Max Exposure is 100 ms

Treatment Plan Field

Description

	Toggle between 180°/360°
---	--------------------------

Buttons



NOTE

- Parameter controls only show on screens when applicable to the pattern selected.
- The Center button requires that the Micro Manipulator is not being used.


Titration

Description

	Select to switch to single spot to perform test burns, adjustable up to 4 spots; select again to revert to the prior laser mode.
--	--

Titrate provides the option to deliver 1 to 4 laser applications in a linear pattern to facilitate the selection of desired laser dosimetry. Power within the titration pattern decreases from left to right, starting with the full power setting and decreasing by one power increment for each spot, moving to the right across the pattern. The titration increment will vary depending on the full power setting and corresponds to the incremental step change in power shown by adjusting the power down arrow; e.g., a 4-spot titration pattern starting at 175 mW will deliver 175 mW>150 mW>140 mW>130 mW.

This feature allows a treatment baseline to be established based on the desired tissue response for the intended treatment and therefore should always be done prior to beginning actual treatment. Initial titration setting levels should be low, progressively moving up to achieve the desired response. Titration burns should only be delivered outside the visual field to the outer retina. It is recommended that the physician wait at least 3 seconds before evaluating the titration lesions because this is the time it takes to achieve the maximum response.

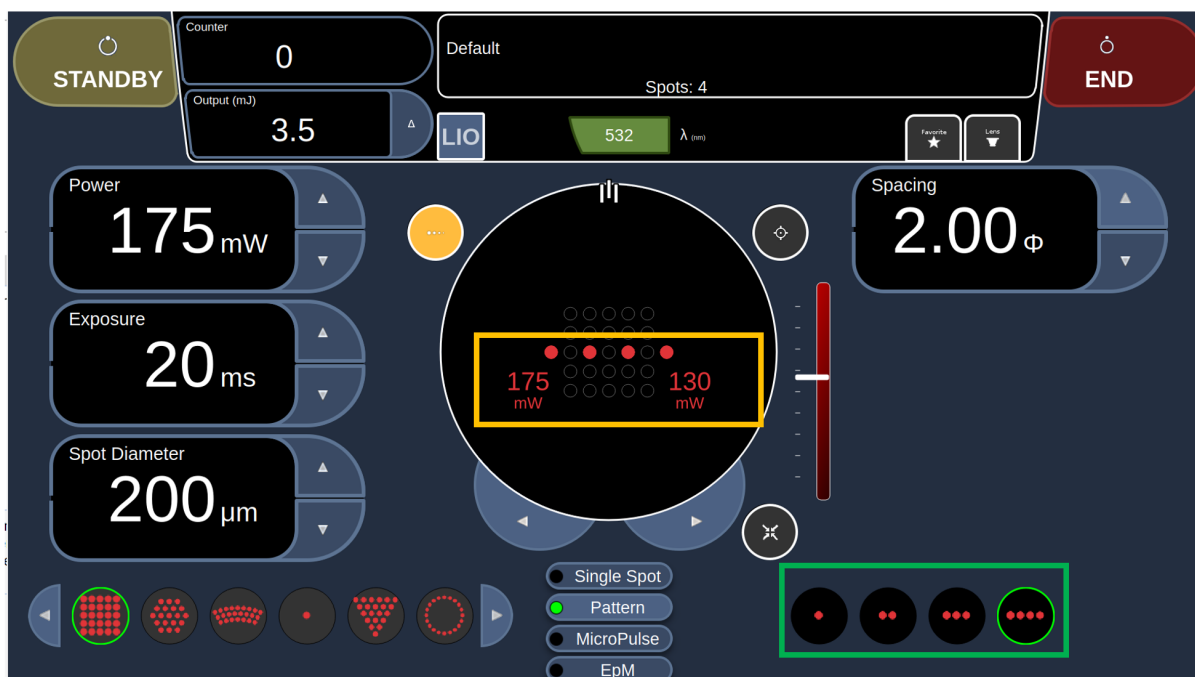

CAUTION

Please confirm that Endpoint Management is automatically disabled and grayed out during titration mode.

Performing titration when Endpoint Management is turned on may result in overpower exposure.



NOTE


The pattern selection bar is disabled when Titration is ON.

Pattern		Starting Point
Pattern	Array 2x2, 3x3, 4x4, 5x5	Single spot at center of pattern
	Rectangle arrays (i.e., 1x5, 2x4, etc.)	Single spot at center of pattern
	Hexagon	Single spot at center of pattern
	Triple Arc	Single spot at center of pattern
	Single Spot	Titration not available
	Wedge 3x3	Single spot at center of pattern
	Wedge 6x6	Single spot at center of pattern
	Arc	Single spot at 12 o'clock of selected ring
	Inner Arc	Single spot at center of selected arc

Pattern		Starting Point
	Triple Ring	Single spot at 12 o'clock in 2 nd ring
	Line (Any # of Spots)	Single spot at center of pattern
	Octant A+B	Single spot at center of selected octants
	Octant A/B (1-7 selected)	Single spot at center of selected octants
	Octant A/B (8 selected)	Single spot at 12 o'clock
	Enhanced Octant Contiguous (1-7 selected)	Single spot at center of the octant to be treated first
MicroPulse: The same as Pattern		
EpM	Enhanced Octant	Single spot at center of the octant to be treated first
	Hexagon	Single spot at center of pattern
	Array	Single spot at center of pattern
ANT	PSLT	Single spot at center of pattern
	Array 3x3	Single spot at center of pattern
	Single Spot	Titration not available

Internal Fixation Light


Description

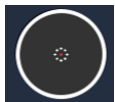	Select to turn on fixation light (when applicable).
---	---

The fixation (guide spot) is available for Octant and Enhanced Octant patterns. The fixation beam is blinking at the center of the aim laser pattern. The fixation beam does not apply to the treatment laser pattern.

Outline


Description

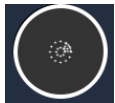	Select to enable full pattern outline (when applicable).
---	--

Display the entire area where the laser irradiation is planned

Auto Advance Feature


Description

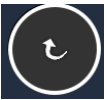	Select to enable automatic advancement (when applicable).
---	---

The selected pattern is delivered when you press the footswitch, after which the system automatically advances to the next pattern.


The orientation is selectable as clockwise or counter-clockwise in Physician Preferences.

System Automatically Advances to Next Pattern

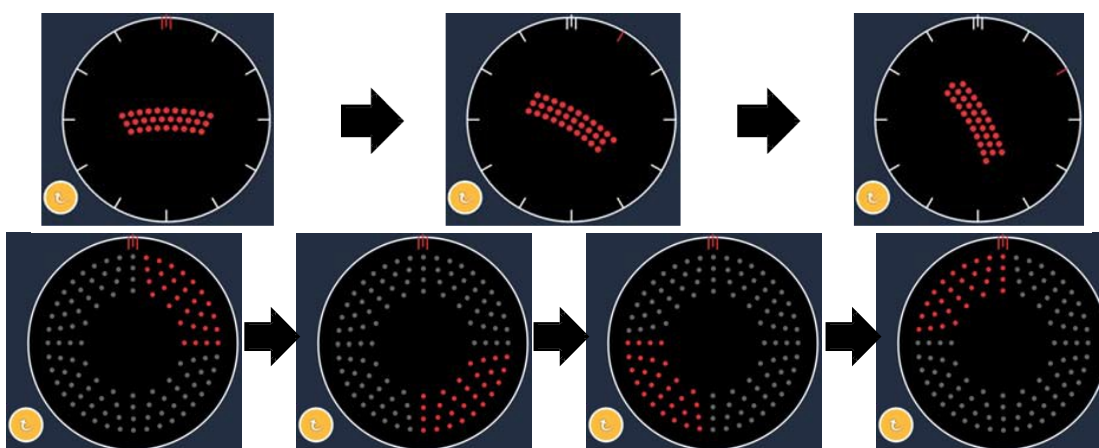

Auto-Transition Feature

If Array or Hexagon pattern is selected in Auto-Transition after Enhanced Octants in the Physician Preferences screen, it allows you to transition to the selected pattern automatically after treating the last subset of EpM Enhanced Octants and to STANDBY mode.

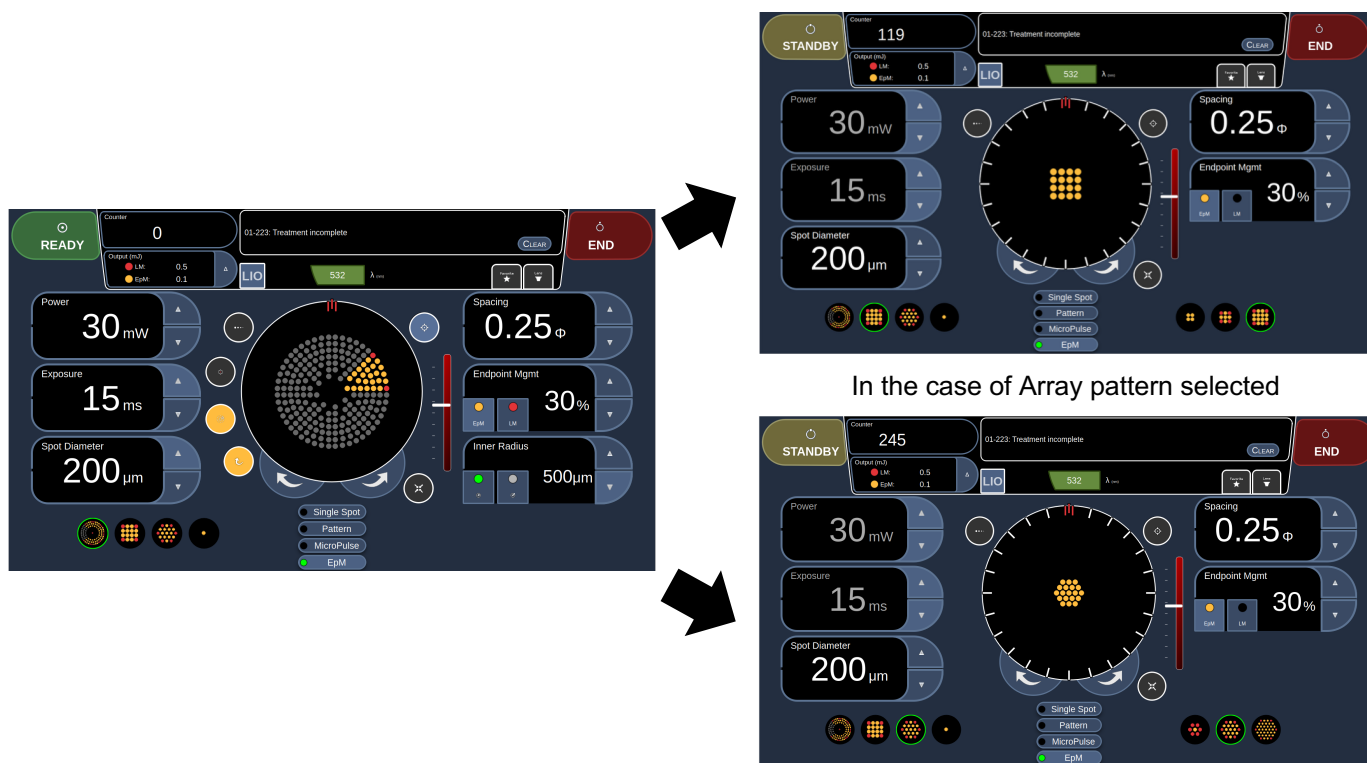


In the case of Array pattern selected

In the case of Hexagon pattern selected

Auto-Transition after Enhanced Octants

AIM

Description

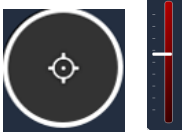	Select to turn on/off aiming beam in STANDBY (always on in READY). Aiming beam intensity: scroll up/down to adjust. Move the white bar to adjust Aim intensity.
---	--

Center

Description


	Select to position the beam to optical center of the treatment.
---	--

Rotation

Description

	Select the arrow to rotate clockwise or counter-clockwise.
---	--

Bottom Panel

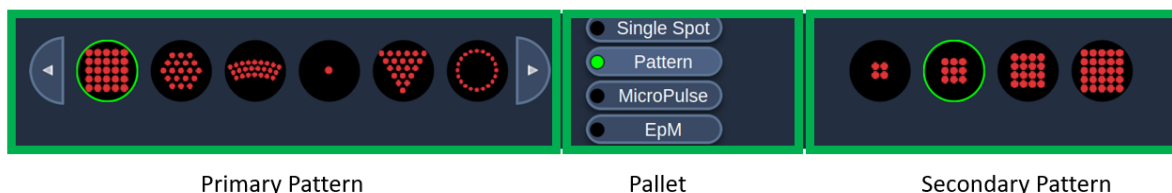

Primary Pattern	Select the pattern type.
Secondary Pattern	Select the pattern size of the primary pattern.
Pattern Pallet	Select the pattern pallet.

Posterior Pattern Descriptions/Pattern Parameters

There are four pallets in Posterior. The choice of pallet determines which parameters (e.g., radius, curvature, repetition rate, etc.) display on the touchscreen control panel, as well as the range of values for each parameter on the Treatment screen.

Pattern Pallet	Primary Pattern
Single Spot	Single Spot
Pattern	Array, Hexagon, Triple Arc, Single Spot, Wedge, Arc, Triple Ring, Line, Octants, Enhanced Octants
MicroPulse	The same as Pattern Pallet
EpM	EpM Enhanced Octants, EpM Array, EpM Hexagon, EpM Single Spot



NOTE

- When you switch pallet among Single Spot, Pattern, MicroPulse and EpM, the setting for each parameter for the pattern will be reset to the default.
- When you select a different pattern within the same pallet, the setting for each parameter will remain the same unless the value is not within the range for the selected pattern. In that case, the setting will be reset to the default.
- The full range of settings is not available for every combination of parameters.

Pattern Pallet	Default Exposure Range	Primary Pattern: Secondary Pattern							
Single	10 ms Default	Single Spot							
Pattern	20 ms Default	Array: 2x2, 3x3, 4x4, 5x5	Hexagon: 7, 19, 39 spots	Triple Arc: N/A	Wedge: 6, 10, 15, 21 spots	Arc: N/A	Triple Ring: N/A	Line: N/A	Single Spot: N/A
	10 ms Default	Octants: A, B, A+B	Enhanced Octants: N/A						
MicroPulse	The same as Pattern Pallet								
EpM	15 ms Default	EpM Enhanced Octants : N/A	EpM Array: 2x2, 3x3, 4x4	EpM Hexagon: 7, 19, 37 spots	Single Spot: N/A				
Anterior	5 ms Default	ANT PSLT: N/A							
	10 ms Default	ANT Single Spot: N/A							
	20 ms Default	ANT Array: N/A							

Single Pallet

Single Spot

Primary Pattern	Secondary Pattern
	Not available


General Use

- PRP, retina tears, retina detachments, iridotomy, trabeculoplasty, LIO

The Single Spot pattern can be used with the slit lamp or an optional laser indirect ophthalmoscope (LIO). Spot diameter and repetition rate settings are adjustable.

Pattern Pallet

Array

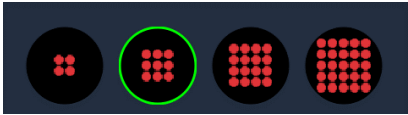Primary Pattern	Secondary Pattern
	

General Use

- PRP, retina tears and retina detachments

The Array pattern is selectable in a variety of shapes and sizes up to a maximum of 25 spots, including rectangle arrays (e.g., 2x3, 4x2, etc.), square arrays (e.g., 2x2, 3x3, etc.), vertical and horizontal lines of up to five spots. Spot diameter and spacing settings are also adjustable.

To select the Array pattern shape and size, drag your fingertip horizontally, vertically, or diagonally across the pattern, or select one of the quick-select buttons at the bottom of the Array pattern screen.



NOTE

Does not allow the user to swipe to single spot for any array patterns. If needed, please select the single spot pattern.

Hexagon

Primary Pattern	Secondary Pattern
	

General Use

- PRP, retina tears and retina detachments

The Hexagon pattern is selectable in three sizes and twenty-four orientations. Spot diameter and spacing settings are also adjustable.

To select the Hexagon pattern orientation, select one of the tickmarks on the rotation guide that displays in the pattern background. The tickmarks are positioned at 15 degree intervals, from 0 to 345 degrees.



NOTE

The selection of spacing will differ depending on the slit lamp used.

Triple Arc

Primary Pattern	Secondary Pattern
	Not available

General Use

- PRP, retina tears, retina detachments


The Triple Arc pattern is selectable in various orientations, depending on the spot diameter, spacing, and curvature settings. You must select the spot diameter, spacing, and curvature settings before selecting the pattern orientation. If you select the pattern orientation and then select the spot diameter, spacing, and/or curvature settings, the pattern returns to its default orientation.

To select the pattern orientation, drag your fingertip along the outside of the rotation guide that displays in the pattern background.

Single Spot

Primary Pattern	Secondary Pattern
	Not available

Wedge

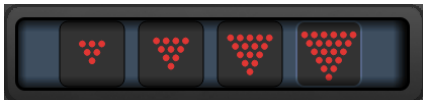Primary Pattern	Secondary Pattern
	

General Use

- PRP, retina tears, retina detachments

The Wedge pattern is selectable in four sizes and twenty-four orientations. Spot diameter and spacing settings are also adjustable.

To select the Wedge pattern size, select one of the quick-select buttons at the bottom of the Wedge pattern screen. To select the pattern orientation, select one of the tickmarks on the rotation guide that displays in the pattern background. The tickmarks are positioned at 15 degree intervals, from 0 to 345 degrees.

Arc

Primary Pattern	Secondary Pattern
	

General Use

- PRP, retina tears, retina detachments

The Arc pattern is selectable in various shapes and sizes, including a single spot, a full circle, and arcs of any number of spots in between. Spot diameter, spacing, and radius settings are also adjustable.

To select the Arc pattern shape and size, drag your fingertip around the Arc pattern, or select the secondary pattern button at the bottom of the Arc pattern screen to select a full circle.

After the outer ring pattern is delivered, the system returns to the inner ring.



System Automatically Advances to Next Ring when Auto Advance Feature is Enabled



NOTE

- If the Auto Advance feature is disabled, only the inner ring displays on the Arc pattern screen.
- The Endpoint Management feature is not available for Arc pattern.

Triple Ring

Primary Pattern	Secondary Pattern
	Not available

General Use

- PRP, retina tears, retina detachments

The Triple Ring pattern is fixed, but the spot diameter, spacing, and radius settings are adjustable. The entire pattern is delivered when you press the footswitch.



NOTE

- The full range of settings is not available for every combination of parameters.
- The Endpoint Management feature is not available for Three Ring pattern.

Line

Primary Pattern	Secondary Pattern
	Not available




General Use

- PRP, retina tears, retina detachments, focal laser

The Line pattern is selectable in nine sizes (2 to 10 spots) and twenty-four orientations. Spot diameter and spacing settings are also adjustable.

To select the Line pattern size, drag your fingertip up or down the Line pattern. Dragging from top to bottom increases the number of spots dragging from bottom to top decreases the number of spots.

To select the pattern orientation, select the **Rotate** button at the bottom.

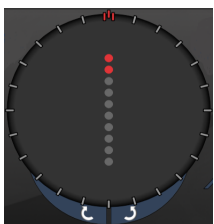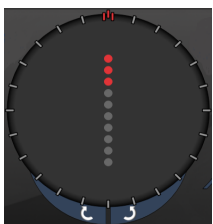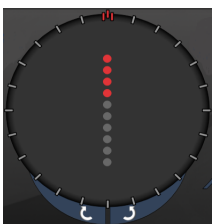
				
1x2	1x3	1x4	1x5	1x6
				
1x7	1x8	1x9	1x10	



NOTE

- The full range of settings is not available for every combination of parameters.
- Do not allow to swipe to single spot for line patterns. If needed, please select the single spot pattern.


Octant

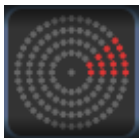Primary Pattern	Secondary Pattern
	


General Use

The Octants pattern is selectable in three formats: **A**, **B**, and **A+B**, as described above. The pattern is selectable in subsets of one to eight octants for A and B formats and one to four octants for the A+B format. Spot diameter, spacing, and radius settings are also adjustable.


To select the Octants pattern format, select the **A**, **B**, or **A+B** button at the secondary of the Octants pattern screen. To select a number of octants, drag your fingertip around the Octants pattern.

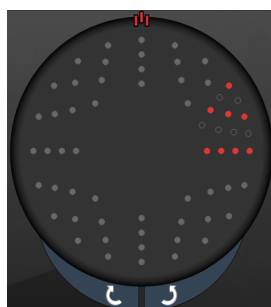

A

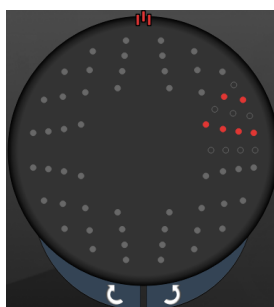

B

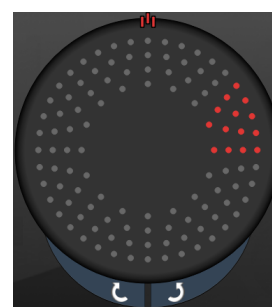

A+B

Enhanced Octant

Primary Pattern	Secondary Pattern
	Not available

General Use

To select a number of octants, drag your fingertip around the Octants pattern.



NOTE

- The selection of inner radius and outer radius will differ depending on the slit lamp used.
- The full range of settings is not available for every combination of parameters.




NOTE

- The LM will be automatically turn off for any auto transitioned pattern after completing the whole EpM enhanced octants pattern. The LM can be enabled if needed.
- When auto-advanced mode is not enabled, there is no auto-transit to configured array or hexagon pattern after completing one segment for EpM enhanced octants pattern.
- When auto-advanced mode is not enabled, the treated segment spots will also be marked as hollow for enhanced octant pattern.



MicroPulse Pattern Pallet

Same pattern selection available as Pattern Pallet in yellow dots.

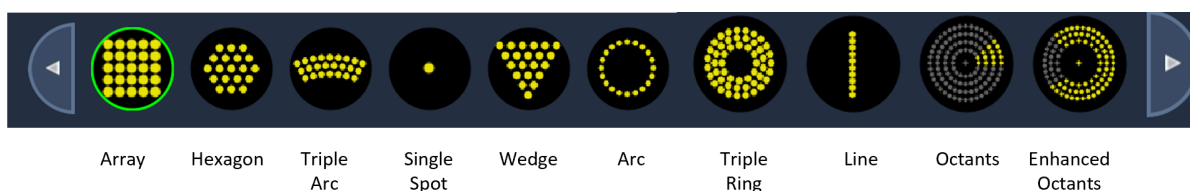

EPM Pattern Pallet

Endpoint Management with Landmark Patterns (optional)

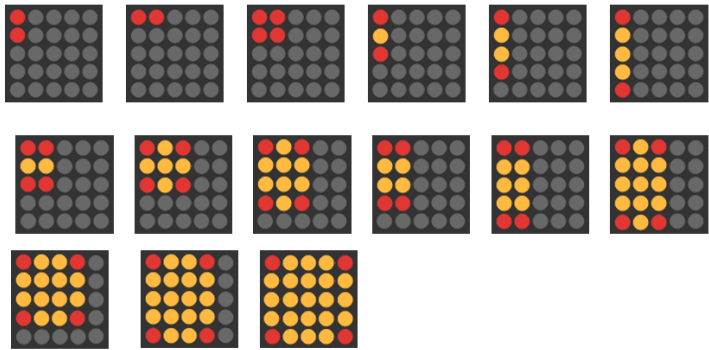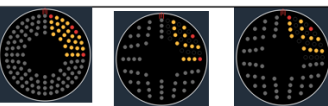Array	
Hexagon	
Triple Arc <Advanced EPM>	
Single	
Wedge <Advanced EPM>	
Line <Advanced EPM>	
Octant <Advanced EPM>	
Enhanced Octant	



NOTE

- The Endpoint Management feature is not available with the LIO.
 - The Endpoint Management feature is only available for EpM ON and LMOFF mode.
-
- The Power and Exposure settings cannot be adjusted for 15ms EpM Single Spot pattern. The Power and Exposure settings keep the same as the previous selected EpM pattern.
 - Endpoint and Landmark patterns only available in Single Spot, Array, Hexagon, Triple Arc, Wedge, Line and Octants pattern pallets when titration is off.


- Establish treatment parameters with Enhanced Octant when in STANDBY mode. The Enhanced Octant pattern will increase the number of “spots” to accompany changes in geometry (spot diameter, spacing, radii). Changes to geometry are only possible when the system is in STANDBY mode. When the system is in TREAT mode, after the first laser application is applied, only changes to Power, Exposure, and Endpoint are allowed.

Anterior Pallet

There are three anterior pattern types.

Pattern Pallet	Primary Pattern
<i>Anterior</i>	PSLT, Array, and Single Spot

PSLT

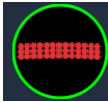
Primary Pattern	Secondary Pattern
	Not available

General Use


- Trabeculoplasty

The PSLT pattern is selectable in various orientations. To select the pattern orientation, select the **Rotate** button at the bottom of the screen. The pattern spacing is fixed, but curvature is adjustable.

Treatment Plan

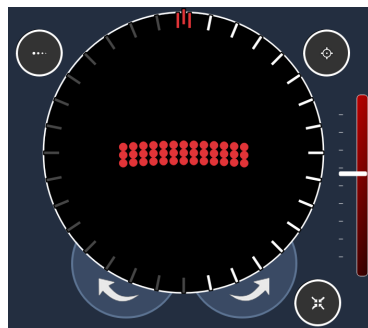


180-degree

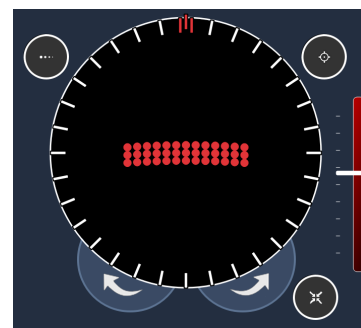

360-degree

Array

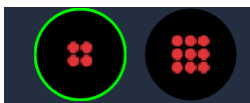Primary Pattern	Secondary Pattern
	

General Use

- Retina tears and retina detachments

The Array pattern is selectable in a variety of shapes and sizes up to a maximum of 9 spots, including square arrays, rectangle arrays, vertical and horizontal lines of up to 3 spots, and a single spot. Spot diameter and spacing settings are also adjustable.

To select the Array pattern shape and size, drag your fingertip horizontally, vertically, or diagonally across the pattern, or select one of the quick-select buttons at the bottom of the secondary pattern screen. To select the pattern orientation, select the **Rotate** button at the bottom of the screen.



NOTE

- *The Endpoint Management feature is not available in the Anterior Treatment.*
- *PSLT is available only for 532 and 577 nm.*

Single Spot

Primary Pattern	Secondary Pattern
	Not available

General Use

- Retina tears and retina detachments

Intraoperative Instructions

Slit Lamp Treatment Procedure

Perform the following procedure:

1. Verify that the slit lamp eyepieces are adjusted to your settings.
2. Position the patient at the slit lamp with their chin on the chin rest and forehead pressed firmly against the head rest.
3. Select the laser treatment spot diameter, exposure time, treatment power level, and pattern type.
4. Position the contact lens on the patient's eye.
5. Select READY mode. The aiming beam will turn on.
6. Adjust the aiming beam intensity.
7. Adjust the pattern spacing, radius, and/or curvature, if applicable.
8. Focus the slit lamp and observe the red aiming beam imaged on the patient's eye. Verify that the laser spots are round and that the pattern is undistorted. Establish proper placement of the laser beam with the slit lamp joystick.
9. Select **Titrate** to switch to a single spot and perform test burns at the periphery of the treatment area.
10. Adjust the laser treatment power for therapeutic effect, and then select **Titrate** again to revert to the selected pattern.
11. Prior to treatment, verify that the power and other parameters are within acceptable ranges.
12. Depress and hold down the footswitch to deliver the treatment laser beam to the tissue.
13. Each footswitch depression will result in **one scanned pattern** unless the treatment is terminated prematurely by early release of the footswitch. Treatment may be interrupted at any point by releasing the footswitch.



NOTE

- *Always place the system in STANDBY mode if there is a prolonged pause in treatment.*
- *If the system is in READY mode and remains idle for 5 minutes, it automatically reverts to STANDBY mode, and the touchscreen control panel becomes light gray. To resume activity, touch the screen. The Status control displays "START UP" for approximately 90 seconds until the system is ready again.*

LIO Treatment Procedure

Perform the following procedure in conjunction with the instructions in the Laser Indirect Ophthalmoscope (LIO) operator manual.

1. Insert adapter into port and verify the LIO software option is activated.
2. Select READY mode. The aiming beam turns on.
3. Select the aiming beam intensity.
4. Select the appropriate exposure time.
5. Select the appropriate laser power.
6. Make any necessary adjustments to the headset to ensure safe and clear view of the retina, as described in the LIO operator manual.
7. Position the aspheric lens and select the spot diameter, as described in the LIO Operator Manual.



NOTE

When using the LIO, the laser spot diameter cannot be adjusted. Spot diameter is determined by the choice of aspheric lens and the position of the LIO relative to the target tissue. See the LIO Operator Manual for detailed information.

8. Depress the footswitch to deliver the treatment laser beam to the tissue.

Between Patient Treatments

At the completion of each patient treatment:

1. Select **End Treatment** to exit the Treatment screen. If the Treatment Report feature is enabled, the Treatment Report window displays; otherwise, the Home screen displays.



NOTE

*If no laser shots were fired during the treatment session, selecting **End Treatment** returns the system to the Home screen whether or not the treatment report feature is enabled.*

2. Disinfect the chin rest and head rest using mild soap and water. Dry with a soft cloth.
3. Disinfect the contact lens following the contact lens manufacturer instructions.

System Shut-down

At the end of the day, or during an extended period of inactivity:

1. Shut down the system as described in "[Shutting Down the System](#)".
2. Remove the key to prevent unauthorized use of the system.
3. Clean the system as described in "[User Maintenance](#)".
4. Place a dust cover over the slit lamp.

Maintenance Instructions

To assure that the system remains safe with regard to electromagnetic disturbance throughout the expected service life, the following actions are recommended:

Annual Maintenance

Preventive maintenance, safety, power and calibration checks is recommended to be performed annually by Iridex Corporation-certified personnel to ensure proper laser performance.

System Repair

All repairs shall be performed by certified personnel to ensure proper system performance.

User Maintenance

The following maintenance procedures are to be performed by the user to ensure proper performance of the system.

Cleaning the Console External Surfaces

Clean the external surfaces of the console daily, after use. Use a cloth dampened with a non-caustic cleaning solution (e.g., soap and water) to clean the external non-optical surfaces of the console. Dry with a clean cloth or allow to air dry. Do not spray or pour cleaning agents directly onto the console.

Cleaning the Control Panel Screen

Use a soft, dry cloth to apply antistatic glass or plastic cleaner to the control panel screen.


Maintaining Effectiveness of Grounds

Clean the unplugged power cord to maintain the protective earth. Use a soft, dry cloth.

Replacing the Fuses

To change the power receptacle fuses:

1. Ensure that the key switch is in the OFF position.
2. Unplug the main power cable from the wall outlet and the system main power receptacle.
3. Insert a small, insulated flathead screwdriver into the fuse holder release receptacle and unlock and remove the fuse holder.

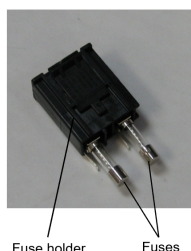

4. Replace the blown fuses with new fuses compatible with the mains supply voltage, as indicated in the [“System Specifications”](#) section of this manual.
5. Replace the fuse holder.

Troubleshooting Guide

If the instrument fails to operate properly, this troubleshooting guide will help you to locate and correct the malfunction. Should a major malfunction occur, contact the Pascal Technical Service Department.

First, check for the following items. If none of these solutions remedies the problem, consult the service manual for additional troubleshooting options.

1. Verify that the wall circuit breaker is in the ON position.
2. Verify that the power cable is correctly attached to the system and the wall outlet.
3. Verify that the key switch is in the ON position.
4. Verify that the door interlock plug is securely connected and, if a door interlock is in use, that the door switch is closed.
5. Verify that the footswitch cable is securely connected.
6. Verify that the LCD control panel cables are correctly installed.
7. Verify that the emergency laser stop button is not depressed.

System does not turn on

Probable Cause:	System is not plugged in.
Suggestion:	Plug in the system. Verify that the power cable is well seated in the wall outlet and in the main power receptacle.
Probable Cause:	Wall circuit breaker is in the OFF position.
Suggestion:	Turn on the wall circuit breaker.
Probable Cause:	The key is absent or in the OFF position.
Suggestion:	Insert key and rotate it to the ON position.
Probable Cause:	Internal system error.
Suggestion:	Turn the key switch to the OFF position, wait at least one minute, and then turn it to the ON position. If the system fails to start, contact Pascal Technical Service.

System monitor is blank for more than 30 seconds

Probable Cause:	Internal system error during boot-up.
Suggestion:	Turn the key switch to the OFF position, wait at least one minute, and then turn to it to the ON position.
Probable Cause:	Monitor/touchscreen panel power, USB and signal cables are not correctly connected.
Suggestion:	Ensure that the monitor cables are correctly connected.

Trouble Adding/Removing Favorites

Probable Cause:	Favorite's Database Update Error.
Suggestion:	Reboot the System.

No aiming beam is present when in READY mode and/or no laser treatment light is delivered when the footswitch is depressed and/or the beams are of poor quality

Probable Cause:	Laser is in STANDBY mode, not READY mode.
Suggestion:	Select READY mode on the control panel.
Probable Cause:	Footswitch is not connected.
Suggestion:	Connect the footswitch.
Probable Cause:	Aiming beam is on low intensity setting.
Suggestion:	Adjust aiming beam intensity on control panel.
Probable Cause:	After five minutes of non-use, system goes to STANDBY.
Suggestion:	Switch mode from STANDBY to READY.
Probable Cause:	Remote interlock has been activated and has disabled system.
Suggestion:	Ensure that the action that activated remote interlock has ceased and proceed.
Probable Cause:	Footswitch and/or footswitch cable damaged.
Suggestion:	Inspect for damage.
Probable Cause:	Internal system error.
Suggestion:	Contact Pascal Technical Service.

Error Messages

When error message appear, select the **Clear** button. If an error code is encountered that is not described within the below table, restart the system (i.e., power cycle the system). If the restart does not resolve the error message, contact Technical Service.

Code	Message	Critical	Clearable	Warning	Description	Action
01	Kill Line	X			The kill lines were asserted in the hardware	Reboot. Call service if error persists.
02	Emergency Stop		X		The emergency stop switch was pressed.	Release Emergency Stop.
03	Footswitch Connect		X		Footswitch has malfunctioned or been disconnected	Connect Footswitch
04	Power Rails		X		An error was detected with the power supplies.	Reboot. Call service if error persists.
05	Watchdog (Galvo)	X			The watchdog fired; this implies the back-end is taking too long in the ISR, or is in an infinite loop.	Reboot. Call service if error persists.
06	Watchdog (Main)	X			The watchdog fired; this implies the back-end is taking too long in the ISR, or is in an infinite loop.	Reboot. Call service if error persists.
07	Interlock		X		The interlock switch was activated.	Connect Remote Interlock
11	Aiming Current Over (High)		X		Aiming beam output above threshold.	Call service if error persists.
15	OPSL Treatment Current Over (Medium)		X		OPSL laser draws electric current above threshold.	Call service if error persists.
16	OPSL Treatment Current Over (High)		X		OPSL laser draws electric current above threshold.	Call service if error persists.
17	Unexpected OPSL Treatment Current		X		Electric current was drawn by OPSL Laser module when it was not expected	Verify device connections.
18	Treatment Current Under for Red 60µm (Low)			X	Red 60µm Laser module draws electric current below threshold.	Call service if error persists.

Code	Message	Critical	Clearable	Warning	Description	Action
19	Treatment Current Over for Red 60µm (Low)			X	Red 60µm Laser module draws electric current above threshold.	Call service if error persists.
20	Treatment Current Over for Red 60µm (Medium)		X		Red 60µm Laser module draws electric current above threshold.	Call service if error persists.
21	Treatment Current Over for Red 60µm (High)		X		Red 60µm Laser module draws electric current above threshold.	Call service if error persists.
22	Unexpected Treatment Current for Red 60µm		X		Electric current was drawn by the Red 60µm Laser module when it was not expected	Verify device connections.
23	Treatment Current Under for Red 200µm (Low)			X	Red 200µm Laser module draws electric current below threshold.	Call service if error persists.
24	Treatment Current Over for Red 200µm (Low)			X	Red 200µm Laser module draws electric current above threshold.	Call service if error persists.
25	Treatment Current Over for Red 200µm (Medium)		X		Red 200µm Laser module draws electric current above threshold.	Call service if error persists.
26	Treatment Current Over for Red 200µm (High)		X		Red 200µm Laser module draws electric current above threshold.	Call service if error persists.
27	Unexpected Treatment Current for Red 200µm		X		Electric current was drawn by the Red 200µm Laser module when it was not expected	Verify device connections.

Code	Message	Critical	Clearable	Warning	Description	Action
28	Base Treatment Temperature Low (Moderate)		X		The base treatment temperature is low	Check fan and ventilation. Allow system to cool.
29	Base Treatment Temperature Low (Severe)	X			The base treatment temperature is severely low	Check fan and ventilation. Allow system to cool.
30	Base Treatment Temperature High (Moderate)		X		The base treatment temperature is high	Check fan and ventilation. Allow system to cool.
31	Base Treatment Temperature High (Severe)	X			The base treatment temperature is severely high	Check fan and ventilation. Allow system to cool.
32	BRF Treatment Temperature Low (Moderate)		X		Laser module treatment temperature is low	Check fan and ventilation. Allow system to warm up.
33	BRF Treatment Temperature Low (Severe)	X			Laser module treatment temperature is severely low	Check fan and ventilation. Allow system to warm up.
34	BRF Treatment Temperature High (Moderate)		X		Laser module treatment temperature is high	Check fan and ventilation. Allow system to cool.
35	BRF Treatment Temperature High (Severe)	X			Laser module treatment temperature is severely high	Check fan and ventilation. Allow system to cool.
36	LBO Treatment Temperature Low (Moderate)		X		Laser module treatment temperature is low	Check fan and ventilation. Allow system to warm up.

Code	Message	Critical	Clearable	Warning	Description	Action
37	LBO Treatment Temperature Low (Severe)	X			Laser module treatment temperature is severely low	Check fan and ventilation. Allow system to warm up.
38	LBO Treatment Temperature High (Moderate)		X		Laser module treatment temperature is high	Check fan and ventilation. Allow system to cool.
39	LBO Treatment Temperature High (Severe)	X			Laser module treatment temperature is severely high	Check fan and ventilation. Allow system to cool.
40	Red Treatment Temperature Low (Moderate)		X		Laser module treatment temperature is low	Check fan and ventilation. Allow system to warm up.
41	Red Treatment Temperature Low (Severe)	X			Laser module treatment temperature is severely low	Check fan and ventilation. Allow system to warm up.
42	Red Treatment Temperature High (Moderate)		X		Laser module treatment temperature is high	Check fan and ventilation. Allow system to cool.
43	Red Treatment Temperature High (Severe)	X			Laser module treatment temperature is severely high	Check fan and ventilation. Allow system to cool.
44	Ambient Temperature Low (Moderate)		X		The ambient temperature is low	Check fan and ventilation. Allow system to warm up.
45	Ambient Temperature Low (Severe)	X			The ambient temperature is severely low	Check fan and ventilation. Allow system to warm up.
46	Ambient Temperature High (Moderate)		X		The ambient temperature is high	Allow system to cool. Reconfigure treatment parameters if necessary
47	Ambient Temperature High (Severe)	X			The ambient temperature is severely high	Allow system to cool. Reconfigure treatment parameters if necessary

Code	Message	Critical	Clearable	Warning	Description	Action
48	Treatment Base Temperature Over Maximum		X		OPSL Base Temperature during treatment exceeds Laser module maximum allowable value	Allow system to cool. Reconfigure treatment parameters if necessary
49	Treatment BRF Temperature Over Maximum		X		OPSL BRF Temperature during treatment exceeds Laser module maximum allowable value	Allow system to cool. Reconfigure treatment parameters if necessary
50	Treatment LBO Temperature Over Maximum		X		OPSL LBO Temperature during treatment exceeds Laser module maximum allowable value	Allow system to cool. Reconfigure treatment parameters if necessary
51	Treatment Temperature Over Maximum (Red Laser modules)		X		Red Laser modules Base Temperature during treatment exceeds Laser module maximum allowable value	Allow system to cool. Reconfigure treatment parameters if necessary
52	Ambient Temperature Over Maximum		X		Ambient Temperature during treatment exceeds maximum allowable value	Allow system to cool. Reconfigure treatment parameters if necessary
56	Aiming Local Light Over (High)		X		Aiming laser output is 100% over the expected value	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service
57	OPSL Treatment Local Light Under (Low)			X	Laser output is 20% under the expected value	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service

Code	Message	Critical	Clearable	Warning	Description	Action
58	OPSL Treatment Local Light Over (Low)			X	Laser output is 20% over the expected value	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service
59	OPSL Treatment Local Light Over (Medium)		X		Laser output is 50% over the expected value	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service
60	OPSL Treatment Local Light Over (High)		X		Laser output is 100% over the expected value	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service
61	Red 60µm Treatment Local Light Under (Low)			X	Laser output is 20% under the expected value	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service
62	Red 60µm Treatment Local Light Over (Low)			X	Laser output is 20% over the expected value	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service

Code	Message	Critical	Clearable	Warning	Description	Action
63	Red 60µm Treatment Local Light Over (Medium)		X		Laser output is 50% over the expected value	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service
64	Red 60µm Treatment Local Light Over (High)		X		Laser output is 100% over the expected value	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service
65	Red 200µm Treatment Local Light Under (Low)			X	Laser output is 20% under the expected value	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service
66	Red 200µm Treatment Local Light Over (Low)			X	Laser output is 20% over the expected value	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service
67	Red 200µm Treatment Local Light Over (Medium)		X		Laser output is 50% over the expected value	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service

Code	Message	Critical	Clearable	Warning	Description	Action
68	Red 200µm Treatment Local Light Over (High)		X		Laser output is 100% over the expected value	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service
69	Unexpected local light (OPSL Aiming)		X		Laser output was detected when it was not expected	Verify device connections.
70	Unexpected local light (Red 60µm)		X		Laser output was detected when it was not expected	Verify device connections.
71	Unexpected local light (Red 200µm)		X		Laser output was detected when it was not expected	Verify device connections.
72	X-Galvo Spot Position Service Error During Aiming		X		An error occurred within the spot position service regarding the X galvo during aiming beam output	Clear error to continue. Call service if error persists.
73	Y-Galvo Spot Position Service Error During Aiming		X		An error occurred within the spot position service regarding the Y galvo during aiming beam output	Clear error to continue. Call service if error persists.
74	X-Galvo Spot Position Service Error During Treatment		X		An error occurred within the spot position service regarding the X galvo during treatment output	Clear error to continue. Call service if error persists.
75	Y-Galvo Spot Position Service Error During Treatment		X		An error occurred within the spot position service regarding the Y galvo during treatment output	Clear error to continue. Call service if error persists.
76	Spot Size Service F-Galvo		X		An error occurred within the spot size service regarding the F galvo	Clear error to continue Call service if error persists.

Code	Message	Critical	Clearable	Warning	Description	Action
77	Software Watchdog	X			Main PCB controller has been reset or did not synchronize with the SBC.	Reboot. Call service if error persists.
78	Host Communication	X			Communication with the SBC failed.	Reboot. Call service if error persists.
79	Port 1 Fault		X		An error was detected with port 1	Reboot. Call service if error persists.
80	Port 2 Fault		X		An error was detected with port 2	Reboot. Call service if error persists.
91	Footswitch Discrepancy		X		Defective Footswitch connected	Verify footswitch connection.
92	Main DSP Communication	X			Communication with the Main DSP failed.	Call service if error persists.
93	Galvo DSP Communication	X			Communication with the Galvo DSP failed.	Call service if error persists.
108	Firmware Revision Request	X			The Main PCB cannot provide the firmware revision	Call service if error persists.
128	Application Configuration	X			The firmware failed to accept the application configuration	Call service if error persists.
129	System Configuration	X			The firmware failed to accept the system configuration	Call service if error persists.
132	Set Audio	X			The firmware failed to enable or disable the audio	Call service if error persists.
204	Touchscreen Error	X			The Touchscreen was disconnected.	Reboot. Check connections. Call service if error persists.
205	Kill Line Asserted	X			The front-end DIO determined that the kill lines were asserted.	Reboot. Call service if error persists.
206	SLA Detached				N/A	N/A
207	Pattern Locked			X	The pattern cannot be modified at this time. Certain patterns can only be modified in Standby mode.	Press Ready button on Touchscreen to return to Standby mode before modifying patterns.

Code	Message	Critical	Clearable	Warning	Description	Action
208	System Busy			X	The user cannot start a treatment because the system has a pending message.	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service
209	Footswitch in Standby			x	The user pressed the footswitch in standby mode.	Press Standby button on Touchscreen to enter Ready mode
210	No output available		X		The user cannot become Ready because the system cannot detect any 1-wire data.	Check ports connection. Reboot. Call service if error persists.
215	Cannot Load Favorite: Not Single Spot			X	The system cannot load the selected favorite because the favorite is not for the single-spot pattern.	Select a different favorite or set parameters using the Touchscreen.
216	Cannot Load Favorite: No Output			X	The system cannot load the favorite because the 1-wire data is not detected.	Select a different favorite or set parameters using the Touchscreen or attach an LIO.
218	Cannot Load Favorite: Bad Pattern			X	The system cannot load the favorite because it refers to a pattern that is not in the system.	Select a different favorite or set parameters using the Touchscreen.
219	Database Error	X			The system cannot function because a database operation returned an error code.	Reboot. Call service if error persists.
220	Laser Warmup Timeout	X			The back-end could not warm the lasers up quickly enough.	Reboot. Call service if error persists.
221	3D Controller Error		X		The 3D controller is disconnected or there was an error attempting to detect it.	Unplug and plug in the 3D controller USB connector to reset. Reboot. Call service if error persists.
222	Invalid Attachment Calibration	X			The calibration file is missing, not readable, or incomplete. Check the log file for the exact problem.	Reboot. Call service if error persists.

Code	Message	Critical	Clearable	Warning	Description	Action
223	Treatment Incomplete		X		This warning occurs when a partial treatment occurs for any reason. (The most common reason is the user released the foot switch too soon, but partial treatments could also be caused by errors.)	User to select End Treatment and Return to the Home Screen. Re-access treatment screen and reconfigure treatment parameters. If error repeats-End Treatment and call Service
225	Cannot Save Favorite		X		Error occurs while saving favorite information to file	Check name provided. Check for duplicates.
226	POST Failed (OPSL)	X			Power On Self-test failure due to OPSL laser initialization error	Reboot. Call service if error persists.
227	POST Failed (Red)	X			Power On Self-test failure due to Red lasers initialization error	Reboot. Call service if error persists.
228	POST Laser Failure				Power On Self-test failure due to both OPSL laser and one or more Red lasers initialization error	Reboot. Call service if error persists.
230	MM controller disconnected	X			Communication with the MM controller failed.	Reboot. Call service if error persists.
231	Invalid 3D Controller	X			An invalid 3D controller is connected.	Power off the system, disconnect the invalid 3D controller from the system and then power on the system.

Calibration Procedure

Regulatory agencies require that manufacturers of US FDA CDRH Class II and IV and European IEC 60825 Class 2 and 4 medical lasers supply customers with power calibration instructions.

Calibration must be performed by an engineer or technician qualified to work on energized electronic laser equipment.

Disclaimer Warning

Calibration of the system is a service procedure to be performed only by certified Iridex personnel. Adjustment by anyone other than certified Iridex personnel voids any existing manufacturer's warranty on the instrument and can result in serious personal injury.

Calibration Instructions

There is no need to open the system for calibration. In fact calibration should be done with the covers on.

Tools Required:

- USB keyboard and mouse
 - NIST-traceable Gentec calibrated optical power meter
 - Power meter to chin rest adapter
1. De-energize the system.
 2. Connect the keyboard and mouse to the USB port.
 3. Attach the power meter to the chin rest adapter.
 4. Install the optical power meter.
 5. Re-energize and start up the system.
 6. Navigate to the Service/Calibration software and enter Service Mode.
 7. Align the optical power meter using the aiming beam.
 8. Calibrate the system as per DC-04077 Service Manual.
 9. Shut down and de-energize the system, as indicated in the "[System Start-up and Shut-down](#)" section in this manual.
 10. Remove the keyboard and mouse.
 11. Re-energize and start up the system.
 12. Verify the calibration using the optical power meter.

If successful, shut down the system and remove the power meter and chin rest adapter. If not successful, repeat the calibration procedure.

System Relocation Instructions

To move the system to another location:

1. Ensure that the wall circuit breaker is turned off and the key switch is in the OFF position.
2. Remove the power cable from the wall outlet and the system main power receptacle.
3. If a remote door interlock is utilized, remove the interlock plug and cable from the interlock port and transport separately.
4. Disconnect the footswitch cable from the footswitch port and transport separately. Never drag the footswitch.
5. Unlock the caster wheel, and position the console a minimum of 15 cm ($\approx$ 6 in) away from walls, furniture or other equipment. Adequate space around the console ensures proper air circulation for system cooling purposes. Lock the caster wheel.
6. If there are changes in environmental conditions (temperature or humidity) allow system to acclimate for 4 hours prior to use.

Room Preparation

- Verify that the system power cable and plug are correctly connected, as instructed in the preoperative instructions.
- Verify environmental conditions are within limits for operation.



NOTE

For systems configured with an electrical wall circuit breaker, always place the electrical wall circuit breaker in the off position before inserting the plug into the receptacle.

- Verify that the electrical service is turned on.
- Verify that the laser warning sign has been posted outside of the treatment room door.
- Ensure that all attending personnel in the treatment room are wearing appropriate eye protection goggles or eyeglasses.

Instructions for Use

Intended Use Environment

The intended Use Environment of the laser system is as follows:

- Professional Healthcare Facility Environment
- Physician offices, clinics, multiple treatment facilities, hospitals except for near active HF surgical equipment and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.

Patient Environment


The Patient Environment is defined by the following:

- When the patient or inspector comes in contact with the devices (including the connecting devices)
- When the patient or inspector touches the person who comes into contact with the devices (including the connecting devices)


In the patient environment (shown below), use the device conforming to IEC 60601-1 standard. If you are compelled to use any device not conforming to the standard, use an insulation transformer that conforms to IEC 60601-1 standard.

Devices applicable to the use in patient environment:

- Insulation transformer

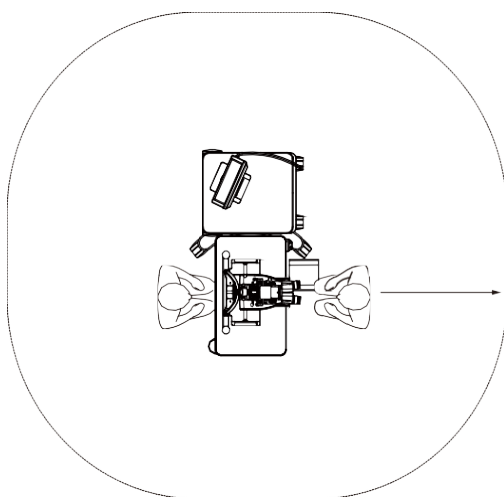

Radius 1.5

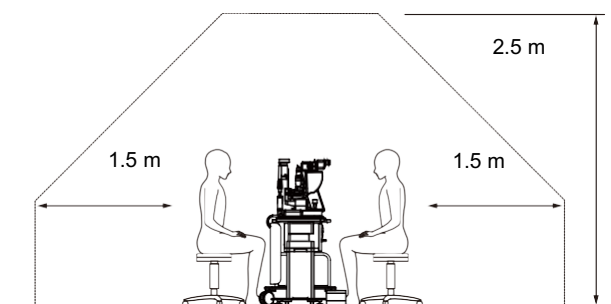

CAUTION

- *Do not use the power strip in the patient environment.*
- *Do not connect an additional power strip or an extension cord to the system.*
- *Do not connect any device which is not recognized as a component of the system.*
- *Use the insulation transformer conforming to IEC 60601-1.*

Electromagnetic Compatibility

Applicable to Iridex Pascal with serial numbers starting with “90” (e.g., 90xxxxxx).

This product conforms to the EMC standard (IEC 60601-1-2: 2014 + A1: 2020). The expected electromagnetic environment for the whole life cycle is professional healthcare facility environment.

- a. MEDICAL ELECTRICAL EQUIPMENT needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in the ACCOMPANYING DOCUMENTS.
- b. Portable and mobile RF communications equipment can affect MEDICAL ELECTRICAL EQUIPMENT.
- c. The EQUIPMENT or SYSTEM should not be used adjacent to or stacked with other equipment. If adjacent or stacked use is necessary, the EQUIPMENT or SYSTEM should be observed to verify normal operation in the configuration in which it will be used.
- d. The use of the ACCESSORY, transducer or cable with EQUIPMENT and SYSTEMS other than those specified may result in increased EMISSION or decreased IMMUNITY of this EQUIPMENT and result in improper operation.
- e. Do not use the devices generating electromagnetic waves within 30 cm from all the parts of the instrument and system. Those devices may have influence on this product.

Item	Cable Shielded	Ferrite Core	Length (m)
ACCESSORIES			
Wireless Footswitch	-	-	-
Laser Indirect Ophthalmoscope (LIO)	-	-	-
CABLES			
AC power cord	Not Used	Not Used	3
USB Cable (for Monitor)	Shielded	Not Used	1
SIP/SOP Cable (for Wired Footswitch)	Shielded	Not Used	4.5

Like other electrical medical devices, Iridex Pascal requires special precautions to ensure electromagnetic compatibility (EMC) with other electrical medical devices. To ensure EMC, your system must be installed and operated according to the EMC information provided in this manual.



NOTE

Iridex Pascal have been designed and tested to comply with IEC 60601-1-2:2014 + A1: 2020 requirements for EMC with other devices.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measure:

- *Reorient or relocate the receiving antenna.*
- *Increase the separation between the equipment and receiver.*
- *Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.*
- *Consult the dealer or an experienced radio/TV technician for help.*

Guidance and manufacturer's declaration - electromagnetic emissions			
The Iridex Pascal is intended for use in the electromagnetic environment specified below. The customer or the user of the Iridex Pascal should assure that it is used in such an environment.			
Emissions test	IEC 60601-1-2 Test Level	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11 A1 Ed. 6.0b: 2016 EN 55011: 2015 + A1: 2017	Class B Group 1 30 – 1000 MHz	Group 1	The Iridex Pascal uses RF energy for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment. The Iridex Pascal is suitable for use in all establishments including domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF emissions CISPR 11 A1 Ed. 6.0b: 2016		Class B	
Harmonic emissions EN 61000-3-2: 2019 IEC 61000-3-2: 2018	Class A Device	per Clause 5 of the standard	
Voltage fluctuations/ flicker emissions EN/IEC 61000-3-3: 2013	Per Clause 5 of the standard	Complies	

Guidance and manufacturer's declaration - electromagnetic immunity

The Iridex Pascal is intended for use in the electromagnetic environment specified below. The customer or the user of the Iridex Pascal should assure that it is used in such an environment.

Immunity test	IEC 60601-1-2 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) EN 61000-4-2: 2009 IEC 61000-4-2: 2008	± 8 kV contact discharge ± 2, 4, 5 & 15 kV air discharge	± 8 kV contact discharge ± 2, 4, 5 & 15 kV air discharge	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst EN/IEC 61000-4-4: 2012	± 2 kV for power supply lines ± 1 kV for input/output lines 5kHz & 100kHz	± 2 kV for power supply lines ± 1 kV for input/output lines Repetition frequency 5kHz & 100kHz	Main power quality should be that of a typical commercial or hospital environment.
Surge Immunity EN/IEC 61000-4-5: 2014	± 1 kV line(s) to line(s) ± 2 kV line(s) to ground	± 1 kV line(s) to line(s) ± 2 kV line(s) to ground	Main power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and Voltage variations EN/IEC 61000-4-11: 2004 + A1: 2017	0% U_T , 0.5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% U_T , 1 cycle 70% U_T , 25/30 cycles 0% U_T , 250/300 cycles	0% U_T , 0.5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% U_T , 1 cycle 70% U_T , 25/30 cycles 0% U_T , 250/300 cycles	Main power quality should be that of a typical commercial or hospital environment. If the user of the Iridex Pascal requires continued operation during main power interruptions, it is recommended that the SL- Iridex Pascal be powered from an uninterruptible power supply or battery.
Power frequency magnetic field EN 61000-4-8: 2010 IEC 61000-4-8: 2009	30 A/m 50Hz/ 60 Hz	30 A/m 50Hz/ 60 Hz	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE: U_T is the a.c. main voltage prior to application of the test level.

Guidance and manufacturer's declaration - electromagnetic immunity			
The Iridex Pascal is intended for use in the electromagnetic environment specified below. The customer or the user of the Iridex Pascal should assure that it is used in such an environment.			
Immunity test	IEC 60601-1-2 test level	Compliance level	Electromagnetic environment - guidance
Conducted Immunity EN 61000-4-6: 2014 IEC 61000-4-6: 2013	0.15 - 80 MHz 3 Vrms & 6Vrms in ISM Band 1kHz AC Mains	0.15 - 80 MHz 3 Vrms & 6Vrms in ISM Band 1kHz AC Mains	Portable and mobile RF communications equipment should be used no closer to any part of the Iridex Pascal, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
Radiated Immunity EN 61000-4-3: 2006 + A1: 2008 + A2: 2010 IEC 61000-4-3: 2006 + A1: 2007 + A2: 2010	10 V/m 80 MHz to 2.7GHz	10 V/m 80MHz to 2.7GHz	Recommended separation distance $d = \frac{6}{E} \sqrt{P}$
Proximity fields from RF Wireless Communications Equipment Immunity EN 61000-4-3: 2006 + A1: 2008 + A2: 2010 IEC 61000-4-3: 2006 + A1: 2007 + A2: 2010	Proximity electromagnetic field from radio communication equipment (See table (a) below)	Proximity electromagnetic field from radio communication equipment (See table (a) below)	where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer, d is the recommended separation distance in meters (m), and E is the radiation electromagnetic field level in volt/meter (V/m).
NOTE 1	These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.		

(a) The table below shows the proximity electromagnetic field from RF wireless communications equipment.

Test Frequency [MHz]	Band ¹ [MHz]	Equipment ¹	Modulation ²	Maximum Output [W]	Distance [m]	Immunity Test Value [V/m]
385	380-390	TETRA 400	Pulse modulation ² 18Hz	1.8	0.3	27
450	430-470	GMRS 460 FRS 460	FM ³ +/-5kHz deviation 1kHz sine	2	0.3	28
710	704-787	LTE Band 13, 17	Pulse modulation ² 217Hz	0.2	0.3	9
745						
780						
810	800-960	GSM 800/900 TETRA 800 iDEN820 CDMA850 LTE Band 5	Pulse modulation ² 18Hz	2	0.3	28
870						
930						
1720	1700-1990	GSM 1800 CDMA 1900 GSM 1900 DECT LTE Band 1,3,4,25 UMTS	Pulse modulation ² 217Hz	2	0.3	28
1845						
1970						
2450	2400-2570	Bluetooth WLAN 802.11 b/g/n RFID 2450 LTE Band7	Pulse modulation ² 217Hz	2	0.3	28
5240	5100-5800	WLAN 802.11 a/n	Pulse modulation ² 217Hz	0.2	0.3	9
5500						
5785						

1) For some devices, only uplink frequencies are included.

2) The carrier shall be modulated using a 50% duty cycle square wave signal.

3) As an alternative to FM modulation, 50% pulse modulation at 18 Hz may be used because while it does not represent actual modulation, it would be worst case.

Ophthalmology References

1. The Diabetic Retinopathy Study Research Group. Photocoagulation Treatment of Proliferative Diabetic Retinopathy Clinical application of diabetic retinopathy-study (DRS) findings, DRS report number 8. Ophthalmology 1981; 88:583-600.
2. Quigley HA. Long Term Follow-up of Laser Iridotomy. Ophthalmology 1981; 88:218- 114
3. Fleischman JA, Swartz MA, Dixon JA. Krypton Laser Endophotocoagulation An Intraoperative Trans-Pars Plana Technique Archives of Ophthalmology, 99: 1610- 1612, 1981
4. Ritch R, Podos SM. Krypton Laser Treatment of Angle Closure Glaucoma Perspectives in Ophthalmology Vol.4, No.4, June, 1980
5. Schwartz AL, et al. Al. Krypton Laser Trabecular Surgery in Uncontrolled Phakic Open Angle Glaucoma Ophthalmology Vol.88 No.3, March 1981.
6. Ophthalmic Laser Therapy Posterior Segment Laser Section Current Techniques Proliferative Diabetic Retinopathy. Edited by Michael J. Bradbury, Vol.1, No.2, 121-129
7. L'Esperance FA, Ophthalmic Lasers Third Edition Photocoagulation of Ocular Disease: Application and Technique. (1988)
8. Kurata F, et al. Intraocular pressure the day of Krypton Laser Trabeculoplasty in Primary Open-angle Glaucoma Ophthalmology 89:338, 1980
9. Thomas JV, Simmons R and Belcher CD III. Complications of Krypton Laser Trabeculoplasty Glaucoma 4:50, 1982
10. Wise JB, Errors in Laser Spot Size in Laser Trabeculoplasty Ophthalmology, 91:186, 1984
11. Kanski JJ, Clinical Ophthalmology, Second Edition. Butterworths, 1989

General information on Usage

Intended Patient Population

The patient who undergoes an examination by this equipment must maintain concentration during the treatment and adhere to the following instructions:

- *Fix the face to the chinrest, forehead rest.*
- *Keep the eye open.*
- *Understand and follow the instructions when undergoing a treatment.*

Intended User Profile

This device must be used by a physician only.

Warranty Information

Iridex Corporation warrants Iridex Pascal to be free from defects in material and workmanship at the original purchaser's location for 12 months.

In order to comply with this warranty, all internal adjustments or modifications must be made by Iridex Corporation-certified personnel or with the express permission of the Iridex Corporation Service Department. The warranty does not apply in the event of misuse, negligence, or accidental damage.

The liability of Iridex Corporation under valid warranty claims is limited to repair or replacement at Iridex Corporation's plant or purchaser's place of business (or, if not practicable, a refund of the purchase price, all at the option of Iridex Corporation).

There are certain other limitations that apply to Iridex Corporation's warranty. Reference should be made to the terms and conditions of sale attached to Iridex Corporation's purchase agreement. Warranty terms may be different by region based on contractual agreements.

Warranty Shipments, Returns, and Adjustments

A warranty claim must be made promptly and must be received during the applicable warranty period by Iridex Corporation. If it becomes necessary to return a product for repair and/or adjustments, authorization from Iridex Corporation must be obtained. Instructions as to how and where products should be shipped will be provided by Iridex Corporation. Any product or component returned for examination and/or warranty repair shall be sent insured and prepaid via the means of transportation specified by Iridex Corporation. Shipping charges for all products or components replaced or repaired under warranty shall be the sole responsibility of the purchaser. In all cases, Iridex Corporation has sole responsibility for determining the cause and nature of failure and Iridex Corporation's determination with regard thereto will be final.

The foregoing Warranty is exclusive and in lieu of all other warranties, whether written, oral or implied, and shall be the purchaser's sole remedy and Iridex Corporation's sole liability on contract or warranty or otherwise for the product. Iridex Corporation disclaims any implied warranty or merchantability or fitness for a particular purpose. In no event shall Iridex Corporation be liable for any incidental or consequential damages arising out of or in connection with the use or performance of the goods delivered hereunder. The essential purpose of this provision is to limit Iridex Corporation's potential liability arising out of this sale.

Disposal

When disposing of the instrument and/or parts, follow local regulations for disposal and recycling or contact your local representative for detailed disposal information.



This symbol is applicable for EU member countries only.
To avoid potential damage to the environment and possibly human health, this instrument should be disposed of (i) for EU member countries - in accordance with WEEE (Directive on Waste Electrical and Electronic Equipment), or (ii) for all other countries, in accordance with local disposal and recycling laws.

This product contains a coin cell.

You cannot replace batteries by yourself. When you need to replace and/or dispose batteries, contact your dealer or Iridex listed on the back cover.



NOTE



EU Battery Directive

This symbol is applicable for EU members states only.

Battery users must not dispose of batteries as unsorted general waste, but treat properly.

If a chemical symbol is printed beneath the symbol shown above, this chemical symbol means that the battery or accumulator contains a heavy metal at a certain concentration.

This will be indicated as follows:

Hg: mercury (0.0005%), Cd: cadmium (0.002%), Pb: lead (0.004%)

These ingredients may be seriously hazardous to human and the global environment.

This product contains a CR Lithium Battery which contains Perchlorate Material-special handling may apply.

See <https://dtsc.ca.gov/perchlorate/>

Decontamination of Returned Equipment

To comply with United States postal and transportation law, equipment shipped to Iridex Corporation for repair or return must be properly decontaminated with a chemical germicide that is commercially available and cleared for sale as a Hospital Disinfectant. To ensure that all equipment has been properly decontaminated, a signed Decontamination Certificate (provided in this section) must be enclosed in the package.

If equipment is received without a Decontamination Certificate, Iridex Corporation will assume the product is contaminated and will assess the customer with decontamination costs.

Any inquiries should be directed to the Iridex Corporation Service Department. These include service of a device, assistance with troubleshooting the device and ordering accessories.

Decontamination Certification

Under the provisions of Postal Law, Title 18, United States Code, Section 1716 and Department of Transportation regulations contained in CFR 49, Part 173.386 and 173.387, "etiologic agents, diagnostic specimens and biological products...are nonmailable..."

The undersigned therefore certifies that the IRIDEX Corporation equipment being returned herein by

Individual/Institution

City, State/Province, Country

Has undergone decontamination with a commercially available germicide cleared for use as a Hospital Disinfectant and is clean and free from biohazards, including – but not limited – human or animal blood, tissue **or** tissue fluids **or** components thereof.

The undersigned also agrees to reimburse Iridex Corporation for any costs incurred in decontaminating the enclosed equipment, in the event said item is received by Iridex Corporation in a contaminated condition.

Model: Iridex Pascal Wavelength: ☐ 532 nm ☐ 577 nm

Serial Number: _____

Iridex Corporation
RMA Number: _____

Position/Title: _____

Name (Printed): _____

Signature

Date (DD/MM/YYYY)

US Technical Service Information

Iridex Corporation
1212 Terra Bella Avenue
Mountain View, California 94043
USA

Phone: +1.650.940.4700

Fax : +1.650.962.0486

PASCALservice@iridex.com